

What's in the Packs

FlowGo Interface Pack



The FlowGo Control Interface



Battery Clip Connector



Power Supply

In Pack
3600 USB-PC



USB Lead with USB Drivers Disk
(suitable for Windows 98SE, Me, 2000 and XP)

OR

In Pack
3600 PC



IBM serial lead

FlowGo Model Pack 1



Fair Model



Model 1 Mimic
Disk



Lighthouse Model



Car Model

Contents

Introducing Computer Control	3 - 7
FlowGo Users Guide	8 - 23
Setting up for the first time and using after the initial set-up	9 - 12
Using Outputs and Inputs	12 - 18
Using FlowGo away from the computer	18 - 19
Power requirements and Technical Specifications	20 - 21
Trouble Shooting	22 - 23
Assembly Instructions for the Models	25 - 30
Installing and Opening the Mimics	31
Introducing the Activities	32 - 36
Keying into the National Curriculum	33 - 36
Teachers Notes and Activity Sheets for the Mimics & Models	37 - 65
Teachers Notes for the Fair Mimic & Model	38 - 43
Fair Activity Sheet	44 - 46
Teachers Notes for the Lighthouse Mimic & Model	48 - 52
Lighthouse Activity Sheet	53 - 54
Teachers Notes for the Car Mimic & Model	56 - 60
Car Activity Sheet	61 - 65

Written in conjunction with:

Mike Martin and Wesley Till

Senior Lecturers at the Faculty of Education University of Central England in Birmingham

Special thanks to:

Flash Ley Community Primary School in Staffordshire

This manual may be copied for use within the premises of the Licensee on condition that it is not loaned, sold or used outside the Licensee's premises.

Copyright: all rights reserved.

Data Harvest's policy is to continually improve products and services, so we reserve the right to make changes without notice. It is acknowledged that there may be error or omissions in this publication for which responsibility cannot be assumed. No liability will be accepted for loss or damage resulting from use of information contained in this manual or from uses as described.

Document No: **DO96**

(Revision 0)

Introducing Computer Control




What is computer control?

Many types of technology can be controlled, including mechanisms and electrical components. Computer control provides unique opportunities for pupils to develop elegant and precise solutions in a simple way. The idea of pupils using a computer to control can appear daunting but in reality when Flowol software is used it is a surprisingly easy rewarding and creative process.

At its simplest, computer control involves using a computer to control electrical devices e.g. lights, buzzers and motors. This may also involve the use of switches as input devices. An interface box, such as FlowGo, is required to interpret the information given to and from the computer and provide the right amount of power for the electrical devices.

It is important to remember that it is we that control the technology, and not the other way around! Most of the technology that we use everyday is simple to control e.g. using knobs and switches to turn things on and off. For those of us involved in educating pupils about control, this is the place to start. The most basic of these involves using a simple switch to control a bulb.

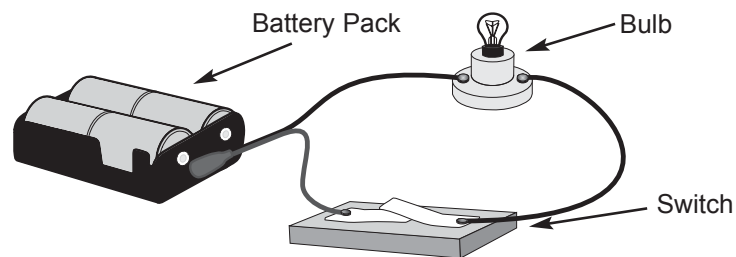

Why use a computer?

Lights, motors and buzzers can easily be controlled by simple on/off switches, magnetic and tilt switches, mechanical timers etc. Why do we need to use a computer?

What makes using a computer better?

Accuracy – especially timing (and positioning)

Capability – controlling more than one thing at once

Power control - motors

Sensing - the environment

Using feedback - so that it is a two-way system

Memory – storing and playing back a sequence of events

Repeating actions – continuously or for a fixed number of times.

Relevance to modern day life

Computer and microprocessor-controlled systems play a significant role in our lives, from traffic management systems and automatic doors to cuddly toys that demand attention. In the modern industrial world, the ability of computers to repeat both easy and complex instructions precisely and endlessly is harnessed to control production lines and to control manufacture to ensure quality products at an economic cost. It is not surprising therefore that increasing attention is being given to computer-aided design, computer-aided manufacture and computer control in design and technology education. A control system is only as good as the person who has designed it! Computer control could be referred to as an invisible technology - you don't notice until it breaks down!



Introducing Computer Control

Computer control and the curriculum

There is long history of using computer control in education. Since the 1960's the programming language LOGO has been used to control 'turtles' on the computer screen and small manoeuvrable robots. Apart from the ability to enable children to move a 'turtle' around, LOGO has a very important role in developing problem solving and thinking skills. Importantly it allows children to 'teach' the computer new things and act upon feedback given from the computer.

Whilst LOGO is relatively easy to use when moving simple 'turtles' on screen it can become complex when controlling a series of motors, lights and buzzers. It has only been since the advent of software that uses flowcharts, such as Flowol, that children of primary age have had real access to computer control.

Currently, computer control is one of the main areas of overlap between Design and Technology (D&T) and Information and Communication Technology (ICT). The present National Curriculum for Design and Technology in the UK includes the use of computer control for children of primary age in both the D&T and ICT Programmes of Study. Now, for the first time, all children are entitled to have experience of using control technology from a young age.

Computer control for a purpose

It is particularly important, with control in design and technology, that a clear purpose for the control activity is established. If the aim of the activity undertaken with children is just to learn how to use the software then it is effectively ICT skills that are being taught. Enabling children to explore a context where a control outcome might be useful will make sure that the software is used as a tool and not just an end in itself. It is where children control a product that they have researched, designed and made that the real links between D&T and ICT are exploited.

What do pupils need to know before starting computer control?

Pupils need two sets of skills:

1. The ability to put together a simple electrical circuit e.g. how to wire up a bulb so they can apply lighting to a model.
2. An understanding of the terminology (sequence of instructions) and the meaning of flowchart symbols. This could be a progression of work started with Learn & Go, LOGO or Roamer. The Flowol mimics could be used to learn how a program is created to perform a sequence of instructions – see Zebra Crossing, and Bridge Lights.

What skills do pupils learn from computer control activities?

Firstly, pupils learn to use the computer for a purpose. At its simplest, pupils learn that computers are not magical and that they need precise instructions in order to perform a task.

Secondly, designing a control system means that pupils are involved in problem solving, and logical thinking.

Lastly, pupils will discover that there can be more than one correct answer
See Key Skills on page 36.

Introducing Computer Control



Computer control in everyday life

Do your pupils recognise control systems when they see them? Ask them to describe what control systems they have used or experienced on the previous weekend. Their experiences may include:

- Traffic lights or pedestrian crossing lights
- Video recorder or interactive T.V.
- Automatic doors at a shopping precinct, supermarket conveyor and bar code reader
- Vending machines for drinks or sweets
- Household appliances e.g. microwave, freezer, fridge, automatic washing machine, or dishwasher.
- Car park barriers, football ground barrier, bus or train ticket machines etc.


Now establish what is common to all of these systems – press a switch, the machine does what it has been 'taught', using a combination of motors and mechanisms, lights and sounds.

Classroom Organisation

Using the I.T. Suite for whole class/large group teaching

Programming skills and a basic understanding of computer control can be taught using Flowol mimics to simulate real-life situations on-screen. Flowol mimics have been designed so a whole class can take part at the same time.

However, a comprehensive understanding of computer control can only really be achieved by using an interface box, such as FlowGo and a model (either ready-made, or preferably built by pupils).

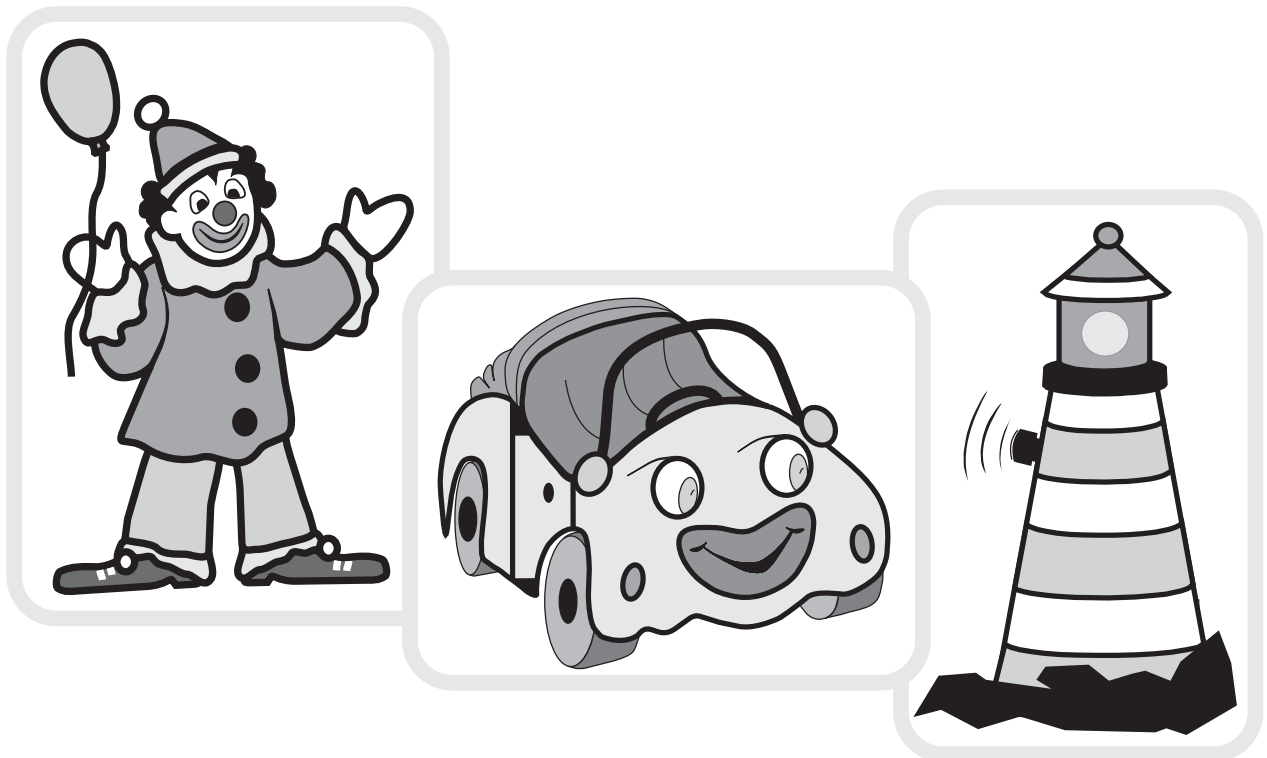




Introducing Computer Control

Stand alone computer

For many teachers, using a single stand-alone machine with an interface attached will be the means of introducing control. In this case more work needs to be done beforehand with children so that their time on the computer is focused. It is important at this stage for the emphasis to be on the story-telling and flowchart creation.

Using models

Ready-made models that match the design of mimics, such as the Lighthouse or Car Alarm are useful in making the link between on-screen simulation and real life control of products. The advantages of using models is that they are instantly available for use and that some guidance on developing the flowchart is provided with them. They can be used for a focused task to develop children's knowledge, skills and understanding in using the control software. It should be stressed that ready-made models are meant for use as an introduction to control. The aim in Design and Technology is for pupil's to develop their own working models. Pupils should have the opportunity to research, design and make their own products using their own choice of mechanisms, electrical components and control systems.

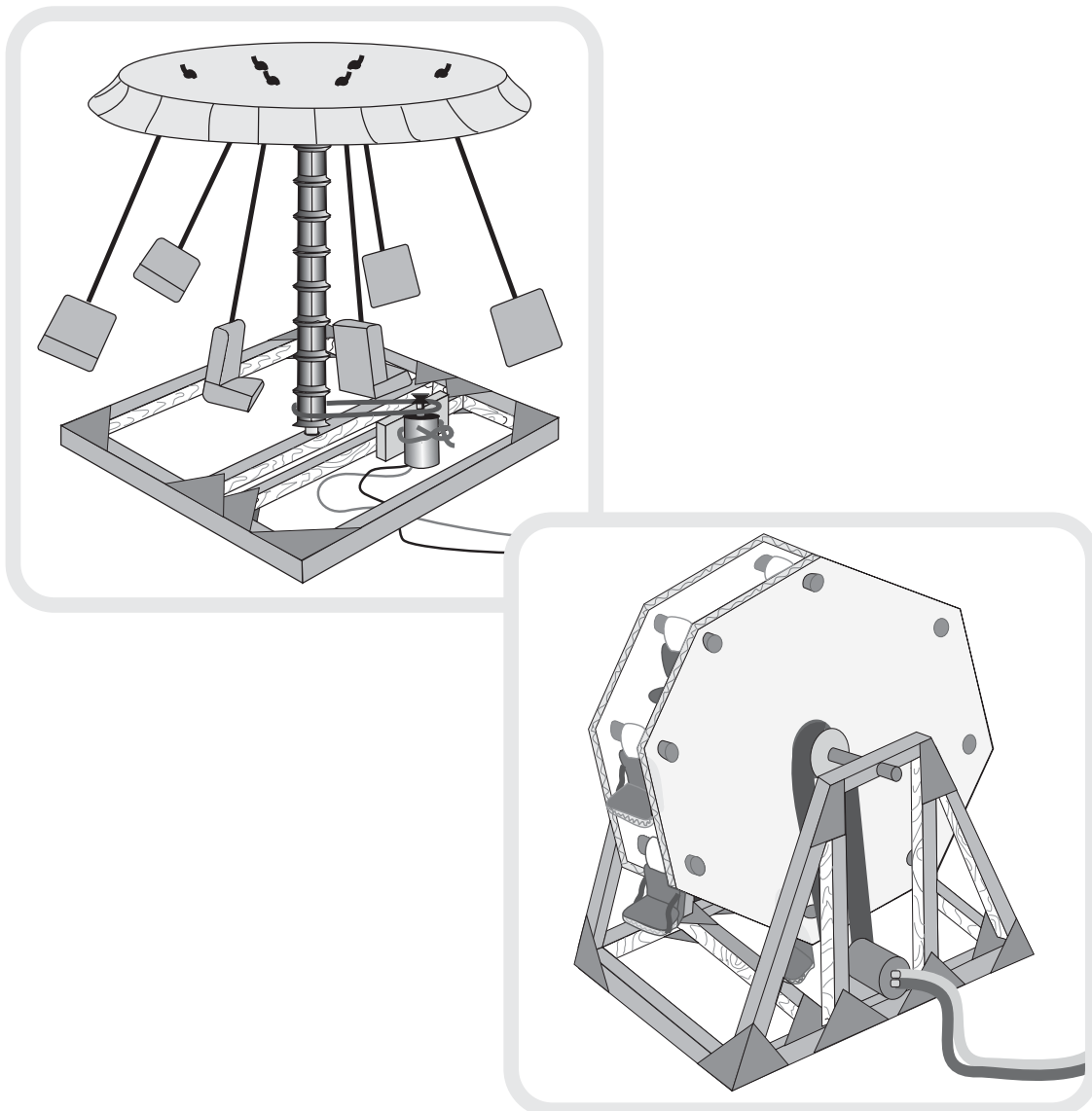

Introducing Computer Control



Moving on – Models with the control system built-in

Have you seen the fridge magnet that meows like a cat, when the door is opened, or the small spinning top that plays a tune whilst it spins.


Pupils at secondary school will be introduced to a whole host of materials to design with; plastics that can be vacuum-formed and cut, metals that can be folded and welded, wood, etc. They will also be introduced to electronics. To make the link with the control that they have used in primary school, Flowol can be used to download a program into a programmable chip, a PIC chip. The information from a control flowchart can be saved onto a chip which is placed on an electronics board and connected to components such as LEDs, motors etc. The pupil's model is then completely independent of the computer. PIC systems enable pupils to design products with a control system that are small and cheap enough to take home!

Conclusion

Simply making a mechanism that can move is an achievement for pupils. Controlling it for a purpose, by using a computer, brings it to life and links what children do in school to their lives outside in a fun, interesting and relevant way.



FlowGo User Guide

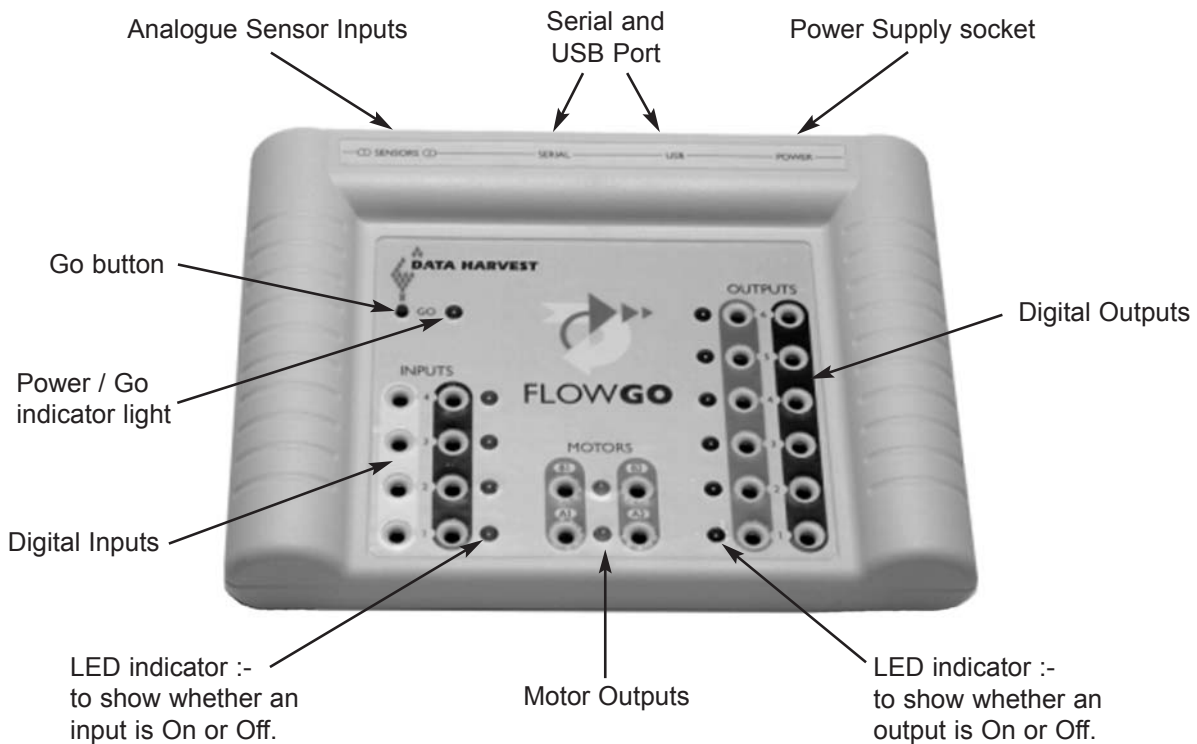

FlowGo is equipped with:

6 Digital Outputs	Red & Black sockets	To power 6-volt electrical devices such as bulbs, buzzers, etc.
2 Motor Outputs	Green sockets	Configured for motors to allow speed control and movement in a forward or reverse position.
4 Digital Inputs	Yellow & black sockets	For switch-type sensors which can register an On or Off state.
2 Analogue Sensor Inputs	Jack sockets at the back of the interface	For level-type sensors capable of detecting changes in environment such as light, temperature, etc.
Power supply socket	At the back of the interface	To power FlowGo either via the mains power supply or the battery snap supplied.
A Serial and a USB port	At the back of the interface	For connecting FlowGo to the computer either via a serial or a USB lead (not both at the same time!).
A green Power / Go light	On the top of the interface	To indicate whether FlowGo is being supplied with power or if a program is running.
A Go push button	On the top of the interface	Used to run the Flowol program stored in FlowGo's memory.



FlowGo User Guide

What is required?

FlowGo is a control interface. An interface is the link between the computer and the electrical components on a model. It provides power for the electrical components and also acts as a buffer between the components and the computer. In order to use a computer to control these components you need Flowol software.

- A computer with a free serial or USB port.
- Flowol software (FlowGo requires Version 2.80 or above).
- The FlowGo control interface.
- A serial or USB connection cable.
- A power source – either use:
 - The mains power supply (supplied with FlowGo) or
 - 4 x size ‘C’ or ‘D’ 1.5v batteries in a 4-cell battery holder and connect via the battery snap connector (supplied).
- A model with electrical components to control.

Setting up for the first time

Turn on your computer.

Step 1

If Flowol is not already on your computer, install the Flowol application. For details of how to install and operate Flowol, please refer to the Flowol tutorial.

Step 2

Connect the power supply lead into the socket on the back of FlowGo and plug into the mains supply. The green LED indicator should stay on.

Note: if the LED flashes, press the Go button once so it stays on.

Step 3


Connect FlowGo to the computer using either the serial cable - see **A** or the USB cable - see **B**.

(A) Using a Serial cable

Insert the serial cable to the serial socket at the back of FlowGo and the other end of the cable to the serial port of your computer.

The serial port on a IBM computer can be either a 9-way or 25-way D-shaped socket and may be labelled as: Serial, COM,RS232

or 

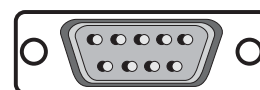

If your computer has more than one serial port, make a note of which port you are connecting to.



FlowGo User Guide

Using a USB cable and its drivers

1. Insert the square-end plug of the USB cable to the USB socket at the rear of FlowGo and the flat USB connector to the USB port on your computer.



Note for Windows 2000 & XP users: the FlowGo drivers need to be installed with administrator privileges.

First Stage:



2. Windows will automatically detect that FlowGo has been attached and will try to locate the drivers for it. A wizard will be generated to help install the drivers for FlowGo. Click on Next.

3. Select 'Search for the best driver for your device' and click on Next.



4. Insert the 'USB Driver For FlowGo' disk (included with the USB cable) into the floppy drive on your computer. Select the floppy disk drive (a:) as the destination for the wizard to search on. Click on Next.

5. Windows will search and then be ready to install the best driver for FlowGo. Click on Next.



6. A progress bar will indicate that the files are being copied.



7. Once installation of the drivers is complete click on the Finish button. Leave the 'USB Driver For FlowGo' disk in the floppy disc drive until the second stage is completed.

Second Stage:

8. The serial drivers that identify FlowGo as a 'virtual' serial device will now need to install. Some operating systems e.g. Windows 98SE, do this automatically (Windows will search, update, and install the FlowGo hardware information). Others systems e.g. Windows 2000, generate a wizard to aid installation: -

- Select 'automatically search for a driver', Next
- Select that the wizard should search on the floppy drive, Next
- The wizard will show it has detected a correct driver from the disc, click Next to start the installation.
- When installation is complete, click Finish to close the wizard.

9. Eject the floppy disc.




FlowGo User Guide

Step 4

To use Flowol with the FlowGo control interface, it must first be set up from the software.

Start the Flowol software application.

1. Leave the **Simulate** option from the Control menu **ticked** until you have made the interface and serial port selection.
2. Select the **Interface** option from the **Control** menu.
3. Select:
 - a. **Serial port**
 - b. the appropriate serial connection e.g. COM1
 - c. the Model of Interface as **FlowGo**
 - d. **OK and Save** for Flowol to save your selection so that it will be automatically configured when next used, or **OK** for these settings to be used for this session only.

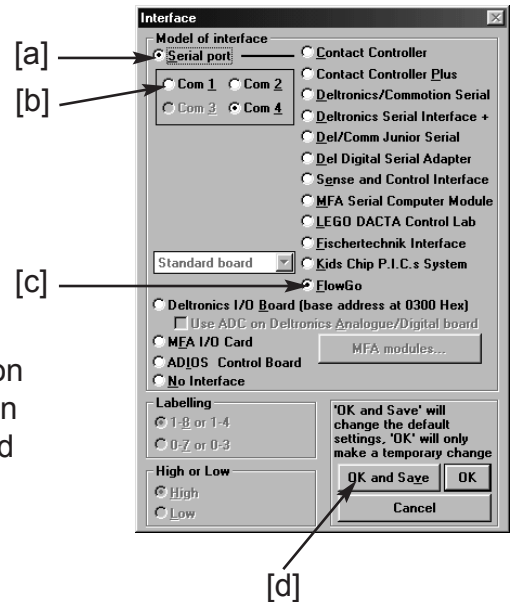

Notes:

- Some ports will be greyed out to indicate that they are not available for selection on your computer.
- If FlowGo is used with a USB lead, it will be seen by Windows as a virtual COM device, so you will need to select the COM port that it has been installed to.
- If you are unsure which number serial port you have used, try each available port in order i.e. Com 1 first, then Com 2, Com 3 & then Com 4. (See Troubleshooting).

A FlowGo button will appear at the bottom of the toolbar on the screen.




Step 5

Untick **Simulate** from the Control menu.

Step 6

To check that communication has been established:-

- Plug an output device e.g. bulb into one of the pairs of output sockets e.g. Output 1
- Select the **Small Monitor** window from the **Window** menu
- Use your mouse to click on the appropriate output number to switch the output device On and Off.

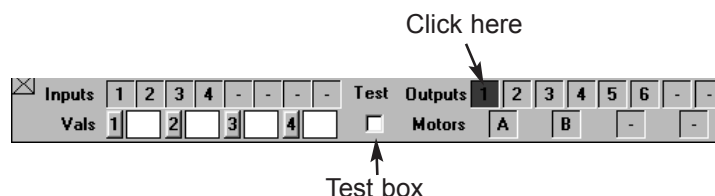

To test an input:-

- Plug a switch to one of the pairs of input sockets e.g. Input 1.
- Click so a tick appears in the **Test box** on the **Small Monitor** window.
- Switch the switch to an On position (e.g. if using a push-to-make switch, press the button) and watch the box below the relevant input light up when the switch is ON.



FlowGo User Guide

Using FlowGo after the initial set-up

- Connect FlowGo to the power supply. The green LED indicator should stay continuously lit.
Note: if the LED is flashing, press the Go button once.
- Connect FlowGo to the computer using the serial or USB cable supplied.
- Start Flowol. If 'OK and Save' was selected during the first set-up, then FlowGo will automatically be configured for FlowGo. If not, re-configure Flowol for use with FlowGo.
- Untick **Simulate** from the **Control** menu.

Simulate

Control programs can be created, tested and modified in Flowol without the control interface attached. This is referred to as 'operating in Simulation mode'.

To simulate FlowGo: -

1. Leave **Simulate** from the Control menu **ticked**.
2. Select **Interface** from the **Control** menu.
3. Select the model of interface as **FlowGo**.
4. Click on **OK and Save**.
5. Flowol will now simulate the correct number of outputs & inputs and the commands available for use with FlowGo.

Using Outputs and Inputs


Connecting Outputs and Digital Inputs:

The Output and Inputs on the front of FlowGo are 4mm sockets.


Methods of connecting electrical components include:

1. Using a pair of 4mm plugs fixed to wire cable, which is then attached to the electrical component (e.g. as used on the components in the Model pack).

Note: These plugs are often referred to as 4mm Banana or Stackable plugs and can have either a solder or screw connection.

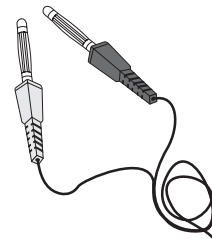

2. Spring Connectors – Product No CK6 – available in packs from Data Harvest.

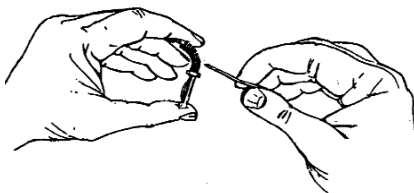

The base of a spring connector is a 4mm plug, the top part is a spring that is pushed back to enable bare wires from electrical components to be inserted.




FlowGo User Guide

Digital Outputs

Plug an output device e.g. bulb, into a pair of output sockets.



Note: If buzzers or LEDs are in use, ensure connection to the correct poles (see below).

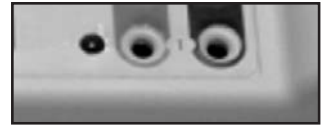

The state of an output socket (On or Off) is shown by the indicator light alongside the socket. When an output is switched **ON** (the line is HIGH) the LED will be **lit**. When the output is switched **OFF** (when the line is LOW) the LED will be **off**.



Examples of output devices include:

Bulbs and bulb holders.

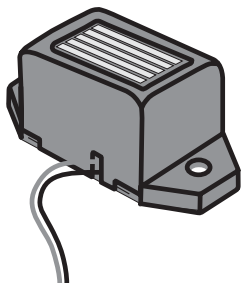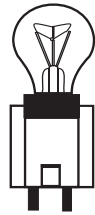

Buzzers:- connect the red lead from the buzzer to the red socket and the black lead to the black socket on FlowGo.

Small D.C. Motors:- can be connected but will only give one direction of rotation.

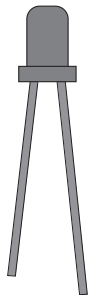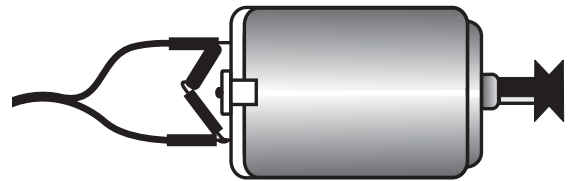

Light Emitting Diodes (LEDs) :- any type can be used. For standard LEDs connect a 470ohm resistor in series with the LED to protect it from the excess current – alternatively a LED with in-built resistor can be used. Connect the long leg of the diode to the red socket and the shorter leg to the black socket on FlowGo.

The voltage from an Output when powered from FlowGo's mains power supply is 6 volts.

The total maximum current supplied by the outputs (6 x outputs + 2 x motors) when used with the mains power supply is 1.3 Amp (1300mA).

Example:

- The current required by 6 x 6v (150mA) ordinary MES bulbs is 900mA, which is within the 1300mA limit.
- The current required by 6 x 6.5v (500mA) High Intensity MES bulbs is 3000mA, which exceeds this limit.

The supply for the digital outputs is via a 'resettable' fuse.

Should the load exceed 1.3 Amps, unplug all the outputs and the fuse will automatically reset after approximately 60 seconds.



FlowGo User Guide


Motor Outputs

The two Motor Outputs can be used with motors to give movement in either direction and/or at varying speeds.

Plug the motor into one of the pairs of motor output sockets i.e. into A1 & A2 or B1 & B2.




Commands in Flowol can be used to change the motors' direction of rotation (to turn it forwards or to reverse it).

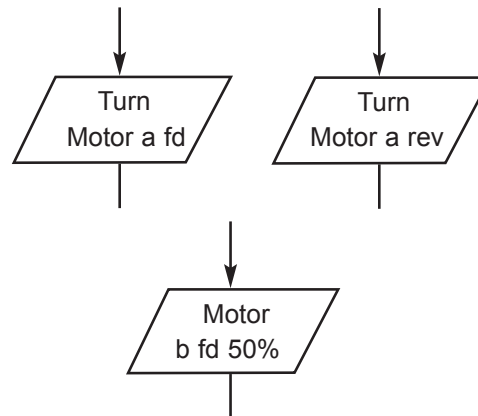

The speed of a motor can be controlled using power commands. If the motor is used without power control, it will run at full speed.

The voltage to the motor output can be pulsed at 10 levels: - 10%, 20%, 30%, 40%, 50%, 60%, 70%, 80%, 90%, 100%.

Select the **Large Monitor** window from the **Window** menu in Flowol. Switch the motor output on and then drag the slider to vary the power to the motor output. This facility is useful when exploring a suitable speed for the motor on a model (providing Flowol is not in simulate mode!).

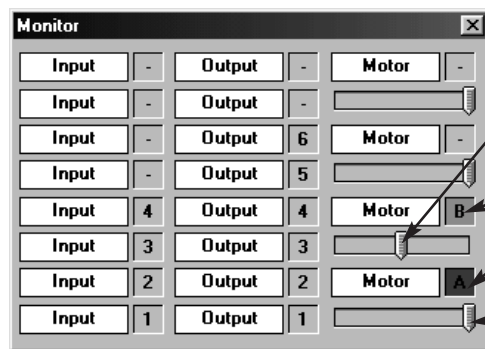


Motor B is turning in a forward direction at 50%

Click here to switch the motor on

Motor A is turning in a reverse direction at 100%

It is possible to use a Motor output as an 'ordinary' digital output whenever extra outputs are required. Use 'Motor' commands to switch the outputs on and off.

The command 'Motor a fd' will switch the left socket (A1) to negative and the right socket (A2) to positive (LED indicator goes green).

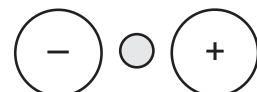

The command 'Motor a rev' will switch the left socket (A1) to positive and the right socket (A2) to negative (LED indicator goes red).




A 'resettable' fuse protects each motor pair. Should the load be exceeded, the fuse will trip. Remove all motor loads and it will automatically reset after approximately 60 seconds.



FlowGo User Guide

Digital Inputs

There are four digital Inputs for use with switches which can register an On or Off state. FlowGo can decide whether or not to carry out a command or run a procedure according to the state of the input.

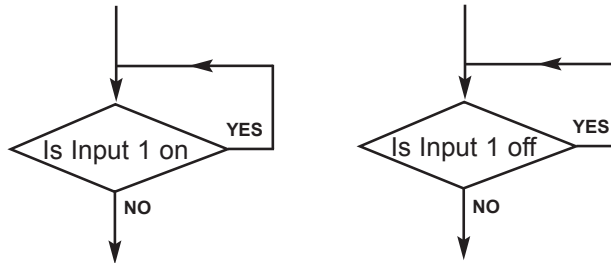


Connect the switch to one of the pairs of input sockets (yellow & black).



Switch the switch to an On position (e.g. if using a push-to-make switch, press the button) and watch the indicator LED next to the relevant input on FlowGo light up when the switch is ON.

A switch is a break in a circuit, which can be joined whenever you wish. Students can make their own switches from everyday materials to suit the own particular needs.

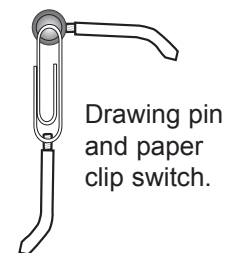


Some examples of the different types of ready-made switches that can be used are:

Push-to-make switch: - when the button is pressed, the contacts connect together allowing the current to flow (Input On): when the button is not pushed, the contacts are open and no current flows (Input Off).



Light Sensitive switch (phototransistor or LDR):- this switch reacts to brightness. If sufficient light is shining at the switch, it passes current and the switch is ON (Input On). When there is not enough light (e.g. by something being placed in front of the LDR), the switch will be OFF. This switch is often used in the construction of a light beam barrier.

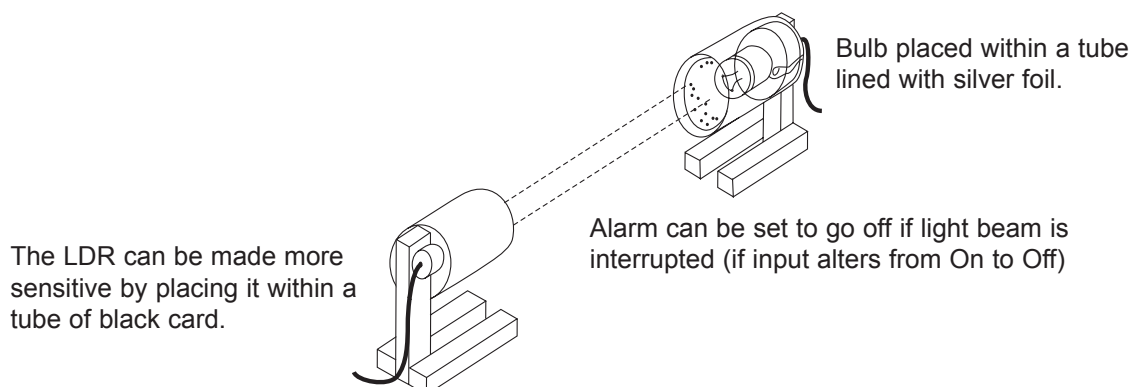




FlowGo User Guide

Tilt switches: - an encased vessel that contains a moving body that makes or breaks a contact when tilted. Tilt the switch up to activated the switch (Input On): tilt it down = no current (Input Off). Often used on a car park barrier.



Magnetic proximity or Reed switch: - used in combination with a magnet. When the magnet gets near, the contacts inside the switch close and the switch is activated (Input On): no magnet = no current (Input Off). Often used to indicate when a door is open or closed.

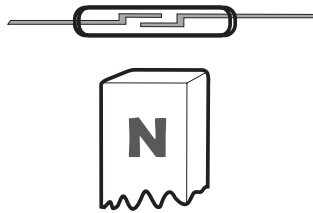


Pressure mat: - a membrane panel switch sealed inside a PVC pad. When there is a load on the pressure mat the switch is activated (Input On): no load = no current (Input Off). Often used in a burglar alarm. Check that the sensitivity of the pad is suitable for the size of the load being used to trigger the switch.




Sensor Inputs

There are two Sensor Inputs for use with analogue sensors e.g. temperature and light level.

Connect an analogue sensor to one of the jack sockets at the back of FlowGo.

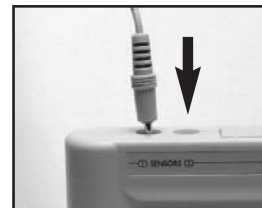

Select either **Large Monitor** from the **Window** menu or **Sensors** from the Graph menu. Click on the down arrow next to the appropriate sensor input number and select the type of sensor from the list.

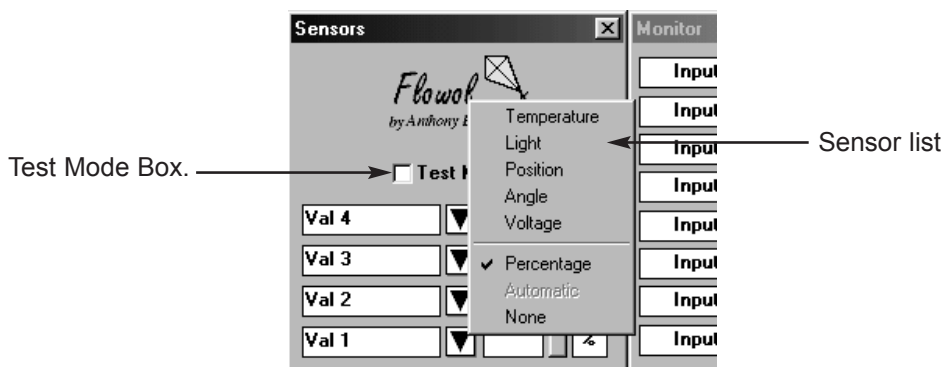




FlowGo User Guide

To see the current value from the sensor, click so a tick is in the Test Mode box (ensure Flowol is not in simulate mode!).

Flowol can decide whether or not to carry out a command or run a procedure according to the value from the sensor.

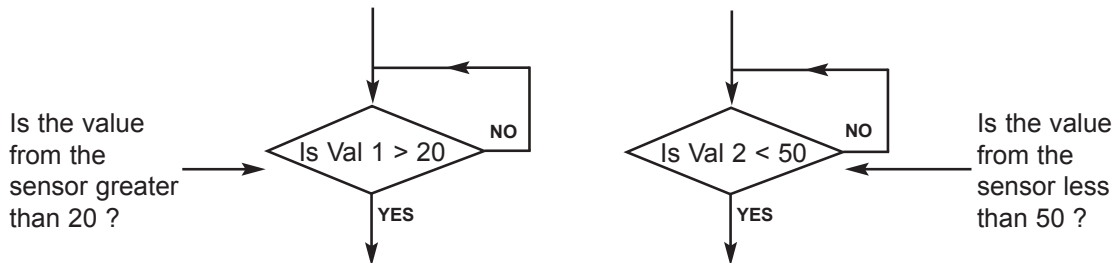


If FlowGo is connected to the computer data is automatically logged whilst a program is running. To view data from the attached sensors, select **Show Graph** from the **Graph** menu.

Click on the appropriate analogue sensor icon in the toolbox. Then click in the Graph window for the graph of that sensor to load.

To load a graph for a second sensor, select its icon from the Toolbox and click on the thick grey line at the bottom of the graph (the gutter).

You do not need to create a flowchart to log data from a sensor, just load the empty graph and click on Run. The timespan defaults to 2 minutes, select Time parameters from the Graph menu to alter. Selecting **Values** allows a line to be scanned across the graph to show analogue values. (See the Flowol Tutorial for more information).

Pre-assembled analogue sensors are available from Data Harvest. If you wish to assemble your own, see 'Technical Specifications'.

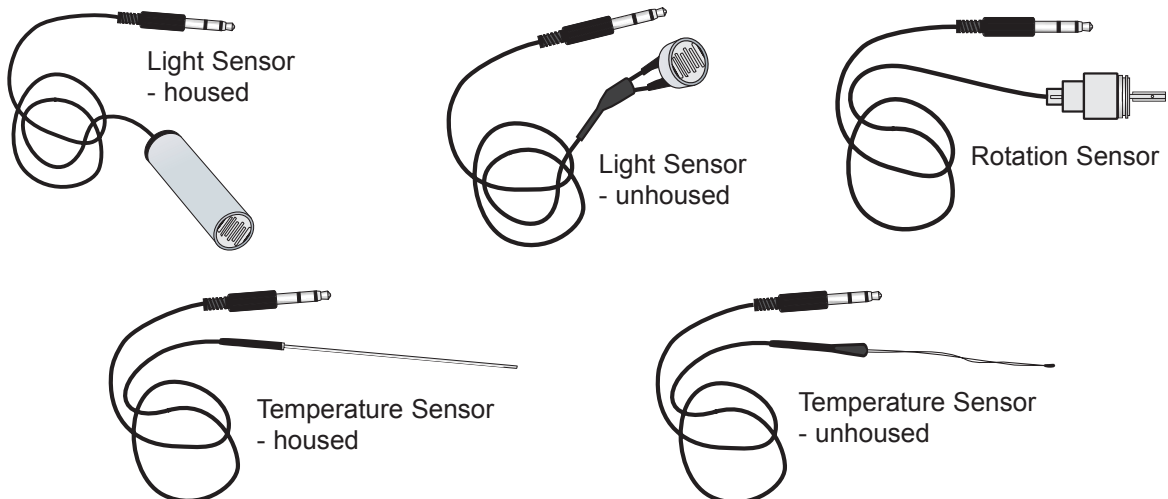



FlowGo User Guide

Using FlowGo away from the computer.

There are two distinct modes of operating FlowGo.

1. It can be used connected to the computer and operated as an ordinary real-time control interface with the Flowol application performing the processing of a program.
2. It can run a program that has been created in Flowol whilst away (disconnected) from the computer. The Flowol program is stored and processed by a microprocessor in FlowGo.

The program is stored in non-volatile memory so even if the power is disconnected or turned off, the program will not be lost. When the power has been reconnected, press Go and the stored program will run.

The microprocessor in FlowGo can be reprogrammed repeatedly.

Downloading a Flowol program to FlowGo

Use Flowol to write your control program in the usual manner - for more details see the Flowol Tutorial supplied with the Flowol software.

Once the control program has been written and tested, it can be downloaded into FlowGo as follows

- Ensure Simulate is unticked and that the control program is not running.
- Check that FlowGo is connected & powered. If the green indicator LED next to the Go button is flashing, press the Go button so it stays on.
- Click on the FlowGo button on the toolbar to initiate the download of the control program into FlowGo.
- As the program is compiled into 'machine code' and downloaded into FlowGo, a progress indicator will be shown.
- Upon completion, disconnect FlowGo from the computer.



Running the Flowol program while disconnected from the computer

Attach a power source (mains adaptor or battery) to the FlowGo interface.

Press the Go button on FlowGo to run the program from start – the green LED will flash while the program is running. When it has finished running, the LED will stay on.

To interrupt the program whilst it is running (green LED flashing), press the Go button once.

FlowGo will hold the compiled program in its memory until a new one is downloaded and it is overwritten.




FlowGo User Guide

Stored program limits

There is room in FlowGo's memory for 256 bytes of program, which is between 100-200 Flowol symbols (the number depends on the complexity of the program). When creating a very complicated program on the computer it is possible to exceed this capacity. For this reason there are limits to the stored program, and they are:-

- **Maximum of 4 Starts.**
If more than four Starts have been used in a program, then only the first four will be downloaded - these will be in order of their addition to the flowchart.
- **A maximum of 16 different Subroutines.**
Example: If (Sub 2) is used five times in a program, then this is counted as one subroutine. If three different subroutines are used e.g. (Sub 1), (Sub 2) and (Sub 3), then these will be counted as three subroutines.

- **Subroutine nesting to 8 levels.** Flowol cannot check for excess levels of subroutine. If the limit is exceeded, the program will appear to download successfully but the program will not execute correctly.

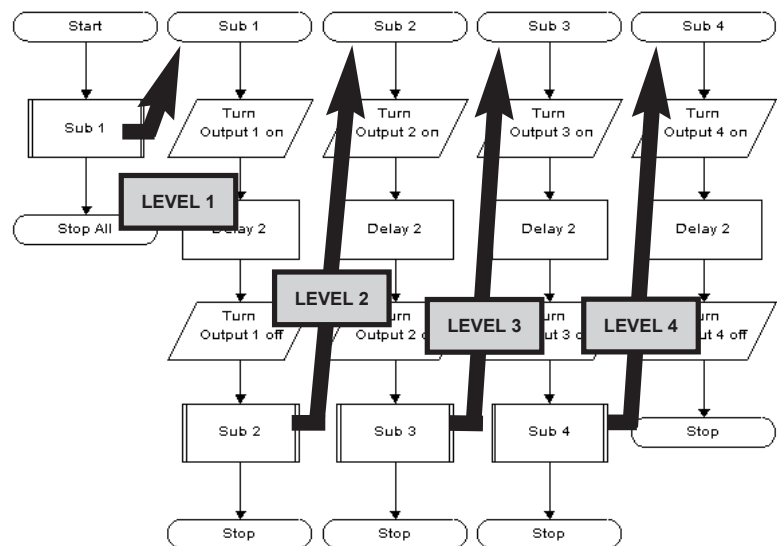

- **No more than 16 different valued Delays.**
Example: If [Delay 2] is used five times in a program, then this is counted as one value of delay. If the delays are for different values e.g. [Delay 1], [Delay 2], and [Delay 3], then this will be counted as three delays.
- **Up to 16 different definitions of a Variable.**
Example: If [Let x = 0] is used five times in a program, then this will be counted as one variable. If [Let x = 0], [Let x = 1], and [Is x = 1] are used then they will be counted as three different definitions.

Note: Only the variables **x** & **y** are available for use with FlowGo.

If the Delay or Variable limits are exceeded, Flowol will not be able to download the program and a 'compile error' message will be displayed.




FlowGo User Guide

Power Requirements

FlowGo has thermal over current protection and will not be damaged by a short circuit if powered using the power supply provided with FlowGo. The thermal over current protection will also prevent excess current being drawn from the power supply or the battery.

Mains Power Supply


Use the power supply packaged with FlowGo to operate it from the mains. It is rated with an output of 1300mA, 6 Volts regulated dc with a positive centre & a negative outer pin.

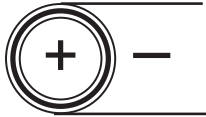

Batteries

FlowGo can be used with a 6-volt battery connected to the power socket using the battery clip connector provided.

We recommend the use of either 4 x size C or 4 x size D 1.5v batteries in a 4-cell holder that has a snap connector. The length of time that FlowGo can operate on a set of batteries will vary according to the age, condition, make of the batteries, the type of components and how much current they use.

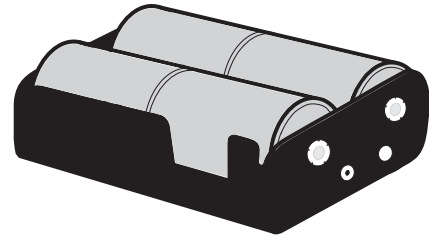

Examples.

1. We used 4 x 1.5v new Duracell Plus Alkaline 'C' batteries to run a program that switched the LEDs on the Fairground entrance model on and off continuously and it ran for 95 hours.
2. We used 4 x 1.5v new Duracell Plus alkaline 'D' batteries to run a program for a buggy with 2 geared down motors going backwards and forwards continuously and it ran for 6 hours.

Safety Notice

Under no circumstances should any part of FlowGo be connected directly or indirectly to any voltage in excess of 6 Volts. Any mains powered supply used to power FlowGo must include a low voltage transformer to fully isolate this unit from the mains supply. The warranty for FlowGo applies only if it has been used with either the power supply provided or with batteries that supply a maximum of 6 volts.



FlowGo User Guide

Technical Specifications

Switching the load to ground drives the digital outputs; each output is capable of handling 100mA, with a combined maximum load from the six outputs of 600mAmps.


The motor outputs can source or handle 200mA continuously. A Pulse Width Modulator, running at a repetition rate of 62 Hz, controls the power.

Digital inputs are pulled high by an internal resistor and configured such that switching an input to ground turns the input 'ON'. The maximum 'ON' resistance is 4K ohms and the minimum 'OFF' resistance is 5K5 ohms. Impedance's between 4Kohms and 5K5 ohms may produce ambiguous results.

The sensors are driven from a single reference voltage of 2.5 V, which can provide up to 1mA of current for each channel.

Sensor pin connections

The connections to the 3.5mm Stereo Jack Socket are as follows:

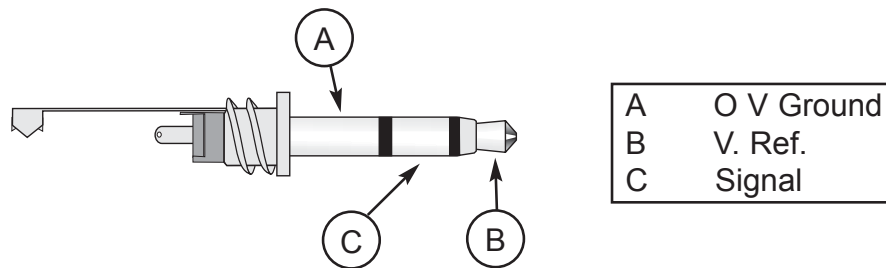

To ensure that 'home-made' sensors respond correctly to the calibration tables found in Flowol, they should be built to the following standards:

Light Sensor

ORP12 (LDR) between VRef. (B) and Signal (C)
2K7 Resistor between Ground (A) and Signal (C)

Rotational Sensor

10K linear pot, track legs to VRef. (B) and Ground (A),
Wiper leg to Signal (C)

Temperature Sensor

100k 6A1 Thermister bead between VRef. (B) and Signal (C)
100k Resistor between Ground (A) and Signal (C)

Environmental Conditions

FlowGo is not waterproof.

FlowGo is suitable for use in a 0 – 30°C operating range and 0 - 95% RH (non-condensing).



FlowGo User Guide

Trouble Shooting

Q. I get an error message ‘ **There is a fault with the serial interface. Reset box and try again**’. This message indicates that Flowol is unable to communicate with the control interface attached.

Note: Flowol does not try to establish communication until you are out of Simulate mode, until you click on the Run button, or until you click on an output or Test an input from the Monitor window.

Check:

1. Is there power to the FlowGo unit? Check to see if the green LED is on.
2. If you are using batteries, disconnect and use FlowGos power supply (the batteries may be drained to such a level that whilst there is sufficient power to light the green power on LED, there is not enough power to allow communication).
3. Is the serial/USB lead connected to the serial/USB port on FlowGo?
4. Is the serial/USB lead connected to the serial/USB port on the computer?
5. Is the LED next to Go flashing – if so, press the Go button so it stays lit.
6. Is the serial port/ Interface selection correct? (See below)

Q. How do I know **which serial port** to select?


1. Make sure Simulate in the Control menu is ticked.
2. Select FlowGo as the Interface from the Control menu.
3. Start with the first Com port that is available (isn't greyed out), e.g. COM1.
4. Click on OK & Save.
5. Deselect Simulate
6. Open the Small Monitor window from the Window menu in Flowol and click on Output 1.
 - If Output 1 on FlowGo lights then your selection is correct.
 - If you get a message to say ‘ **There is a fault with the serial interface. Reset box and try again**’, then repeat the above from Step 1 – 6 but select the next Com port that is available. (Make sure Simulate in the Control menu is ticked before you begin the selection process).

Q. I have switched the **Output on** and the device attached is **not working**: -

1. Does the screen display indicate that the output is switched on?
2. Is the red LED next to the Output lit?
3. Are you still in Simulate mode?
4. Is FlowGo attached to a power source?
5. Is there a loose connection?
6. Is the bulb loose in the bulb holder?
7. Have you connected the leads from your buzzer or LED correctly?

Q. When I select Interface from the Control menu **FlowGo is not in the list**.

1. Select About from the Help menu. The version number of Flowol is on the bottom right side of the flash screen. The version number needs to be 2.80 or higher.
2. If you have an earlier version of Flowol, contact Data Harvest for information about an upgrade.

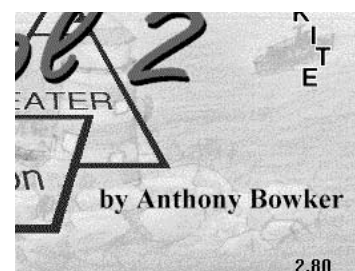



FlowGo User Guide

Q. I have the Flowol **upgrade**, but I already have Flowol on my computer. Should I delete the old version first?

1. Yes. Exit Flowol if it is open.
2. Select Uninstall from the Start, Programs, Flowol 2 menu.



Q. When I open a Monitor window or a Mimic, I don't get the **right number of outputs and inputs** available.

1. Go to Interface from the Control menu and check that the interface selected is FlowGo.

Q. When I try to download to FlowGo, I get an error message 'Failed to connect to FlowGo – check connection and try again (**Error No 300**)'.

1. Is the LED next to Go flashing – if so, press the Go button so it stays lit.
2. Is FlowGo powered up, and the serial lead connected?

Q. When I try to download to FlowGo **nothing happens**.


1. Check that Simulate is unticked.

Q. When I run my control program I get a message to say '**Error parsing decision**'.

1. You are downloading a program that has been created using variables other than those available for FlowGo e. g. (n)
2. Alter the program using either x or y as the variable.
3. Alter the Interface selection to FlowGo so that in future the options available reflect only those available for FlowGo.

Q. How can I check which serial port the USB driver has loaded to?

1. From the Start menu, select Settings, Control Panel, and System Properties.
2. Select the Device Manager tab
3. Select "Ports (COM & LPT).
It will be shown in the list e.g. FlowGo [COM 4].

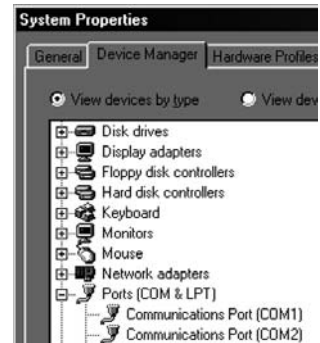

Q. How do I uninstall the USB drivers?


1. The USB lead should not be connected.
2. From the Start menu, select Settings, Control Panel, and Add/Remove programs
3. Select the Install/Uninstall tab
4. Select 'FlowGo USB Drivers' from the list.
5. Click on the Add/Remove button.

Warranty

FlowGo is warranted to be free from defects in materials and workmanship for a period of 12 months from the date of purchase provided it has been used in accordance with any instructions, under normal conditions. This warranty does not apply if the FlowGo has been damaged by accident or misuse.

In the event of a fault developing within the 12-month period, FlowGo must be returned to Data Harvest for repair or replacement at no expense to the user other than postal charges.

Assembly Instructions for the Models

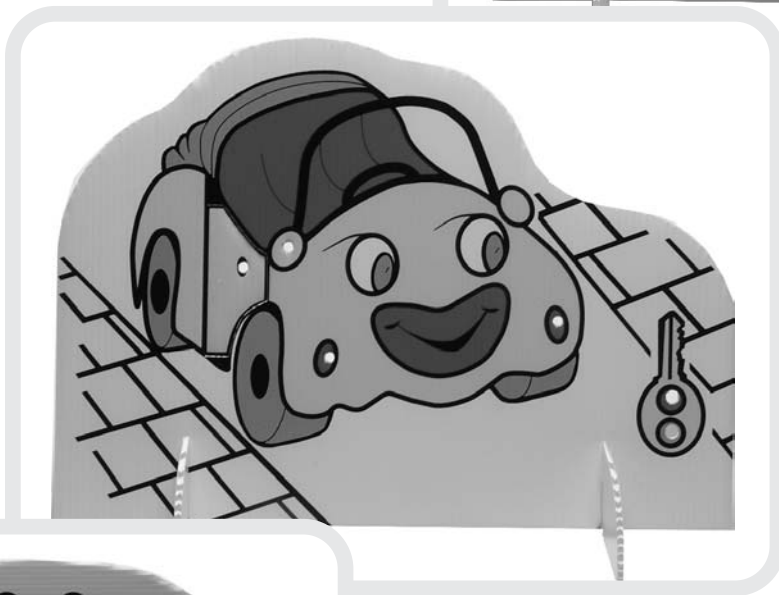


Assembly Instructions for the Fair Model




The Model comes with all the electrical components pre-wired with colour-coded 4mm plugs (red and black) to match the Outputs on FlowGo.

The contents of the pack are:-

- A 12-way ribbon cable with 6 coloured light emitting diodes (LEDS),
- One double-sided sticky pad.
- The model (with its stand) showing a clown at the entrance to the fair.

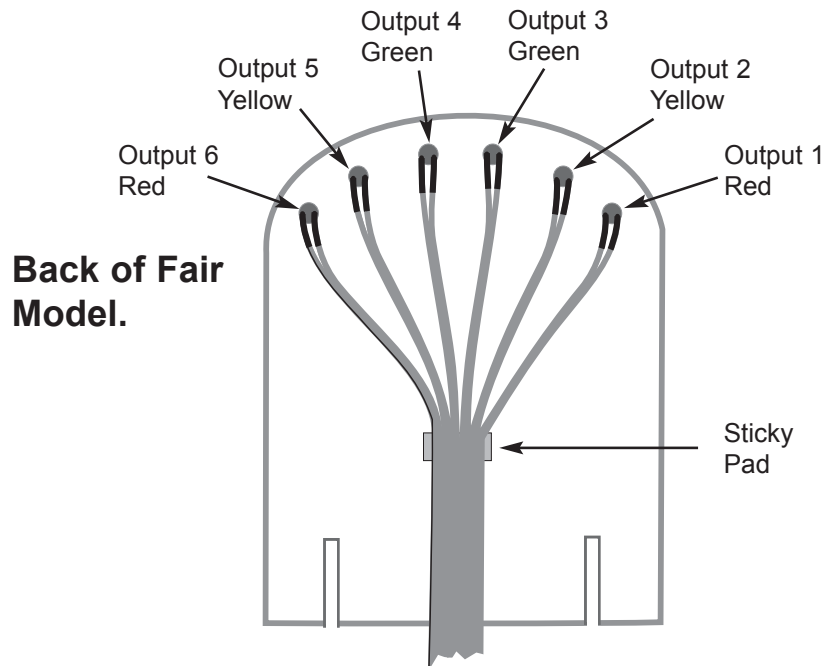

1. The LEDs need to be fed into the holes from the back of the picture. Hold the LED securely just below the diode and push gently into the holes (take care not to bend the lead/diode junction).

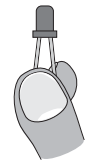

2. You may find it useful to label each pair of output plugs to ensure the LEDs 'match' the labelling for the mimics. Work along the cable from the back of the model from left to right.

The red LED is output 1.
The yellow LED is output 2.
The green LED is output 3.
The green LED is output 4.
The yellow LED is output 5.
The red LED is output 6.



3. Use the double sticky pad (provided) to fix the ribbon cable to the back of the model.


4. Assemble the stand (see page 30)




Assembly Instructions for the Lighthouse Model

The model comes with all the electrical components pre-wired with colour-coded 4mm plugs that match the Output / Inputs on FlowGo - red and black for Outputs, yellow and black for Inputs.


The contents of the pack are:-

- A 6-way ribbon cable that has one lamp, one buzzer, and a connector,
- A light dependent resistor (LDR) this has a clear lens through which it is possible to see an orange snake-like shape. 
- Two double-sided sticky pads.
- The model (with its stand) showing a lighthouse.

Assembly

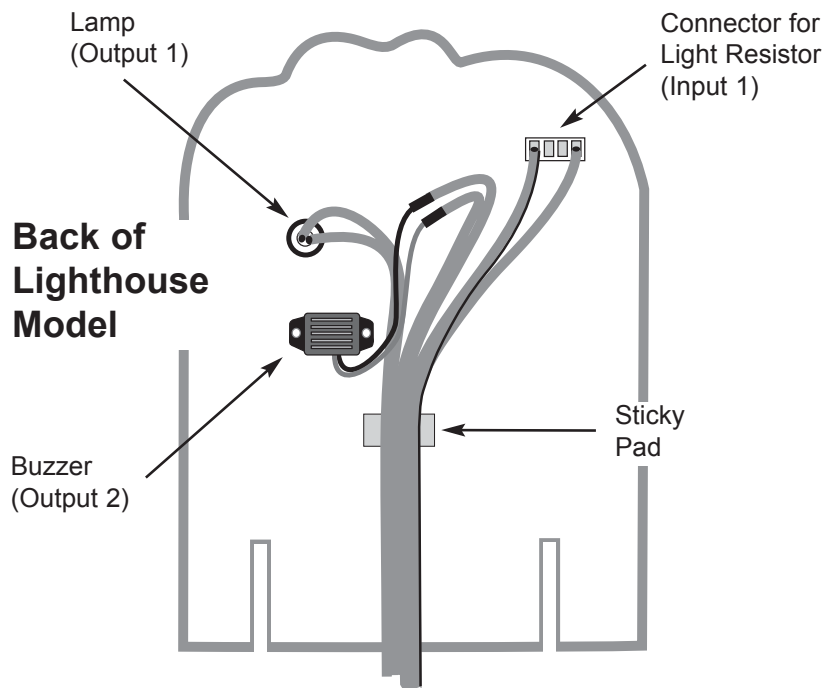


1. Push the lamp into the hole from the back of the model until its holder just makes contact with the front of the model (take care not to bend the lead junction).
2. Use a double sticky pad to secure the buzzer into position on the rear of the model.
3. Feed the thin metal legs of the light resistor from the front of the model through the slot cut in the Sun (see picture **A**).

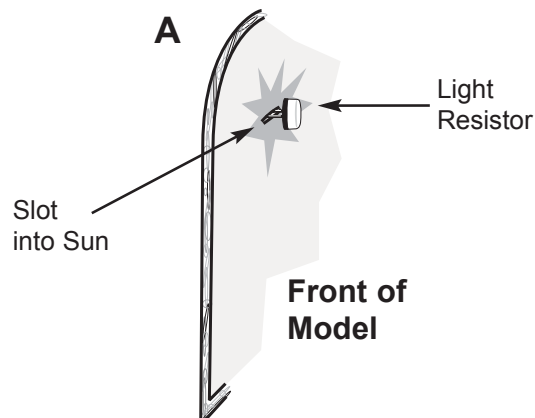

Assembly Instructions for the Lighthouse Model




4. Insert the legs from the resistor into the two outer holes on the connector (picture **B**)

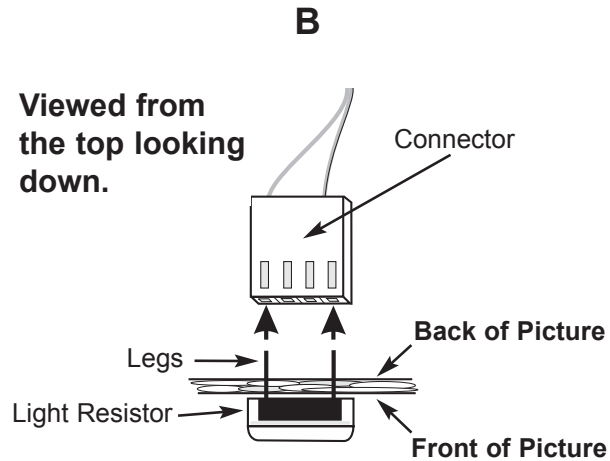


5. You may find it useful to label each pair of plugs to ensure they 'match' the labelling on the mimics. Work along the cable from left to right from the back of the model;

The buzzer is output 2.
The lamp is output 1.
The light resistor is input 1.

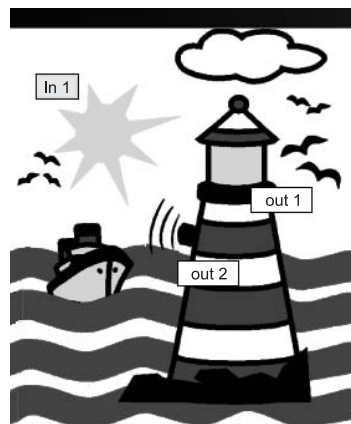

6. Use a sticky pad to fix the ribbon cable to the back of the model.
7. Assemble the stand (see page 30)

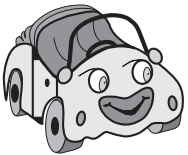

Assembly Instructions for the Car Model

The model comes with all the electrical components pre-wired with colour coded 4mm plugs that match the Output / Inputs on FlowGo - red and black for Outputs, yellow and black for Inputs.

The contents of the pack are:-

- A 14-way ribbon cable that has one red LED, two yellow LEDs, one magnetic proximity switch, (with magnet), one buzzer and two push-to-make switches,
- Four double-sided sticky pads,
- The model (with its stand) showing a car .
- Rivet and knob for door handle.


Note: You will need scissors to cut the sticky pad and may find a small spanner or pliers useful for the push switch assembly.

Assembly


1. Push the yellow LEDs into the correct holes from the back of the model.
2. Use a double sticky pad to secure the buzzer onto the rear of the picture roughly in the position shown.
3. Unscrew the small hexagonal nut and washer from one of the push switches. Push the switch into the correct hole from the back of the model. Screw the nut and washer into place from the front of the model to secure.

Note: You should be able to get a sufficiently tight fit using your fingers, if not use a small spanner or pliers.


4. Repeat with the other push switch.

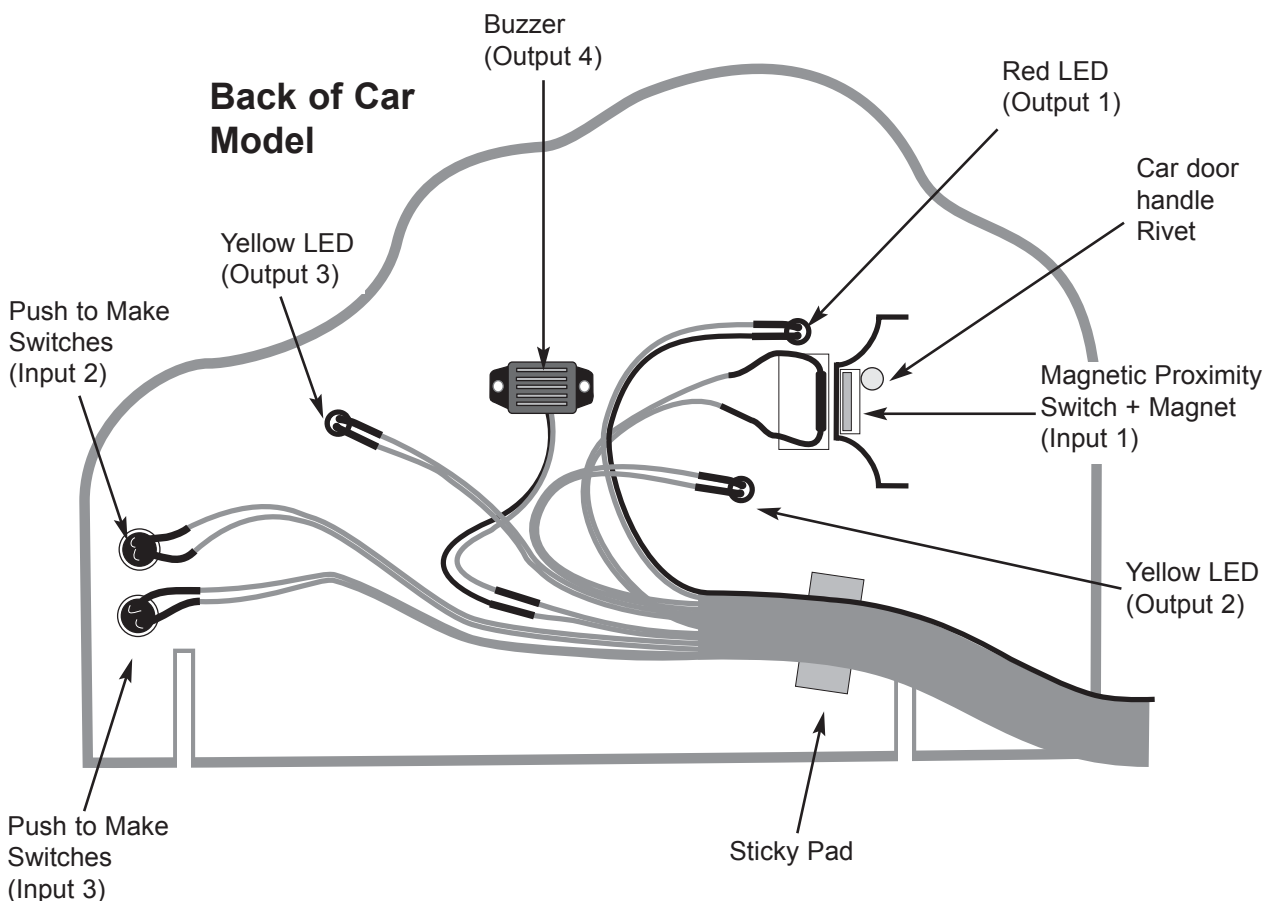


Assembly Instructions for the Car Model

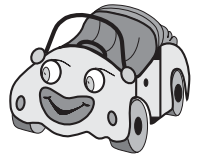

5. Cut one of the double-sided sticky pads to a suitable size for the magnet (20mm x 6mm). Stick the pad **(1)** close to the edge of the door on the back of the model as shown on sketch **A**. Stick the magnet onto pad **(1)**.
6. Stick the other sticky pad **(2)** on the edge of the door frame – keep the pad below the hole for the red LED and close to the edge. Press the magnetic proximity switch on to pad **(2)**. Do not bend the lead/switch junction.
 Note: it is necessary for the magnet and switch to be close together when the door is closed for the switch to operate correctly.
7. Push the red LED into the correct hole from the back of the model.

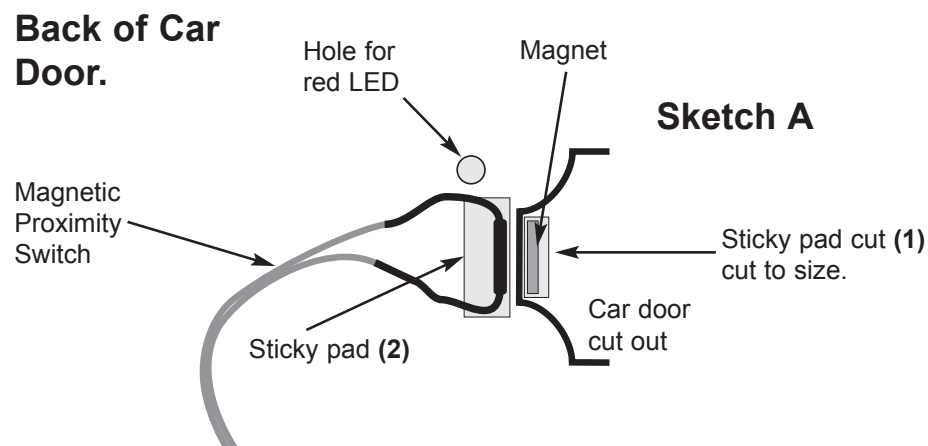


8. Insert the plastic rivet into the hole for the door handle. Hold the rivet in position and fix the handle to the door from the front of the model (as shown in sketch **B**)

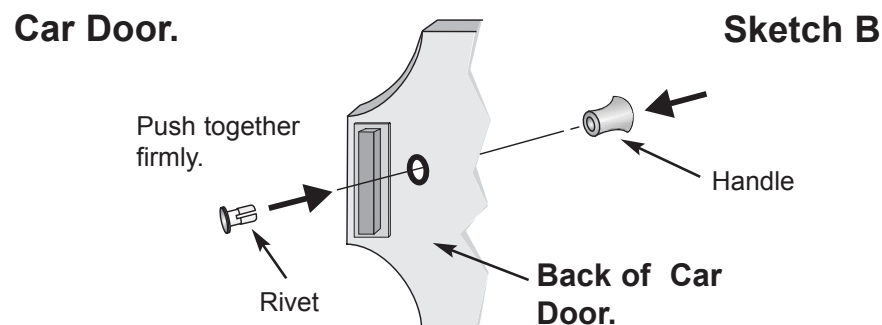

9. You may find it useful to label each pair of plugs to ensure they 'match' the labelling for the mimics. Work along the cable from the back of the picture from left to right;


- The first push switch is input 3.
- The second push switch is input 2.
- The buzzer is output 4.
- The magnetic proximity switch is input 1.
- The yellow LED is output 3.
- The next yellow LED is output 2.
- The red LED is output 1.

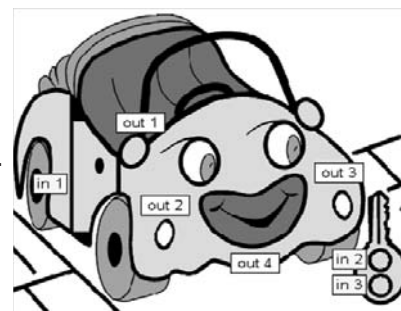


10. Use a sticky pad to fix the ribbon cable to the back of the model.
11. Assemble the stand (see page 30)



Assembling the Stand

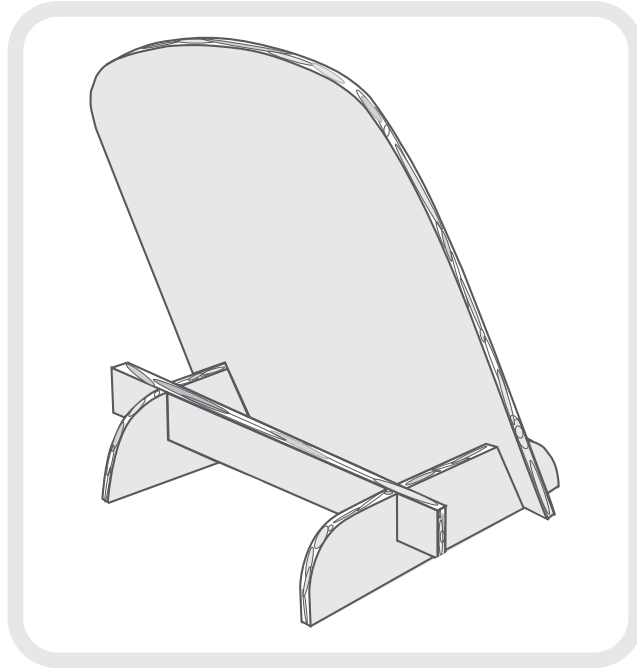

- Assemble the stand and rear support to the main model slots.

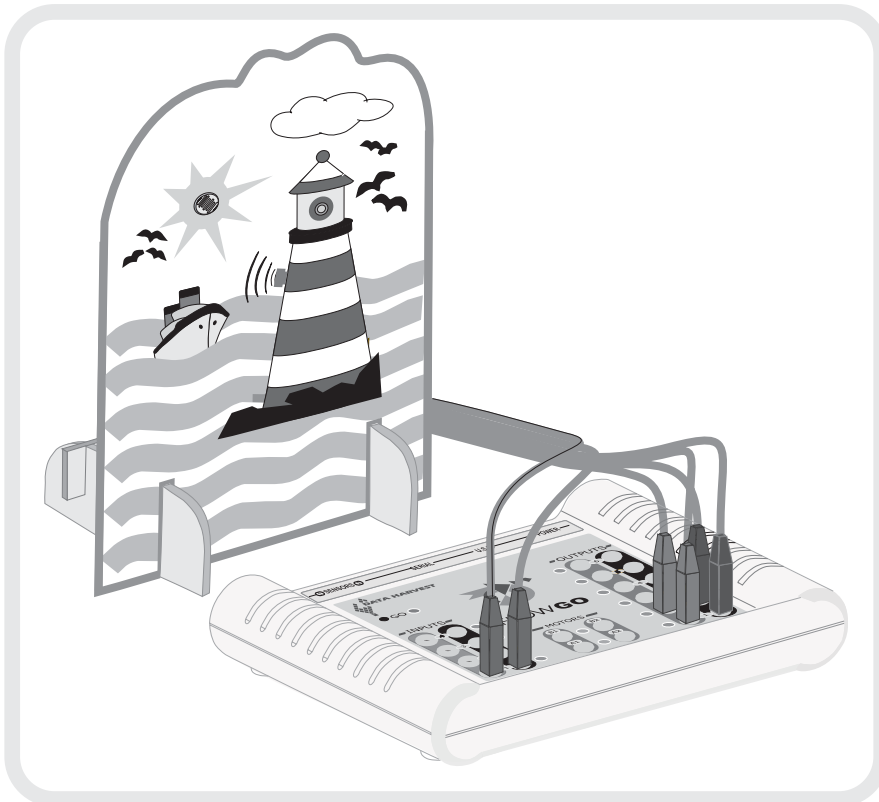

- You are now ready to plug the wires into the Outputs and or Inputs on FlowGo.

Installing and Opening the Mimics



Installing the Model Pack 1 Mimics.

Note: Flowol2 must be installed on the computer before the mimics can be installed.

PC Windows:

1. Place the floppy disc labelled 'Model Pack 1 Mimics' in drive A.
2. From the **Start** menu select **Run**.
(Win 3.1 Users select Run from Program Manager)
3. Type **a:setup** in the dialogue box and click on OK.
4. Follow the on-screen instructions. The installation utility will automatically search for the mimic folder that was created when Flowol 2 was installed on your hard disc i.e. 'Mimic folder found at C:\DHG\Flowol\Mimics'. Click on OK to install the new mimics.



Note: If a copy of the Flowol program was open (in use) before installation started, exit & restart before trying to open the mimic.

5. When installation has finished eject the Mimic disc.


Opening a Mimic

Start Flowol.

Click on the **Window** menu and select **Mimic**. Select the appropriate mimic from the list.

The names for the mimics in this pack have been prefixed by p1.

Tick in the Show labels box to view the input/outputs labels.

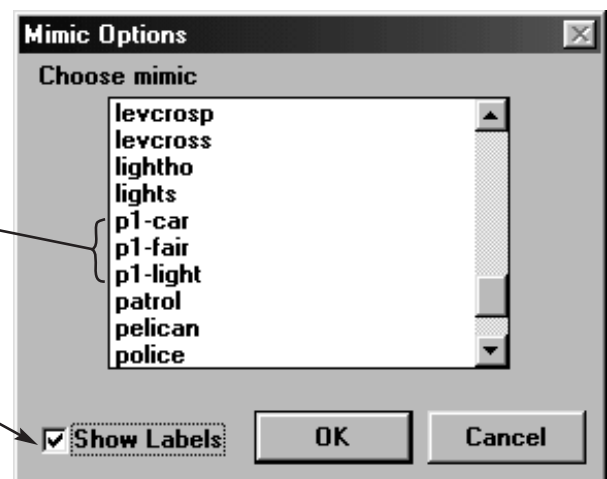



Introducing the Activities

The following teacher's notes and activities are designed to introduce computer control with Flowol software and the starter models in Model Pack 1. They aim to give a progression through the following skills:

- Creating a program that turns a single output on and off.
- Creating a program that turns more than one output on and off.
- Creating a program that runs continuously
- Using an input to start a program.
- Using a sub-routine.
- Repeating a sub-routine for a set number of times.
- Using an input to call up a sub-routine.
- Using more than one flowchart.
- Using more than one input.
- Using inputs to stop, as well as start a sequence of actions.
- How to tackle more complex programs.

These notes and activities have been written to suit primary teachers with pupils aged 9 – 11, as starter activities for computer control. It is assumed that teachers will have read the Flowol Tutorial and made themselves familiar with the Flowol program before starting these activities.

The Teacher's Notes give an overview of the pupil's activity sheet, with teaching points for each activity, plus example programs for all the challenges. There are suggestions for extension activities that refer to other Flowol mimics available separately. These activities have been included to help those with mixed ability pupils, who wish to set further challenges to the achievers or skills practice. A list of the additional Flowol mimics are on the following page.


Each Activity Sheet shows several activities, each with an example flowchart, plus challenges for pupils to work through independently, or in pairs or groups. Each challenge is designed to use the skills taught in the previous activity.

Introducing the Activities



Additional Flowol Mimics

Primary Mimics	Robot Face, Pelican Crossing, simple Railway Level Crossing and simple Automatic Home
More Primary mimics	School Patrol Crossing, Police Car, Fire Engine, Water Chute, Kitchen and simple Home Alarm System
Secondary mimics	Pelican Crossing, Railway Level Crossing, Car Park Barrier, Greenhouse and Automatic Home
More Secondary mimics	Water Chute (with counter), Kitchen, Home Alarm (with keypad), Theme Park Water Flume and Jam Factory
Control Pictures	Traffic Lights, Burglar Alarm and Washing Machine
Commotion Mimic Disc	Clown Face, Giant Traffic Lights, House

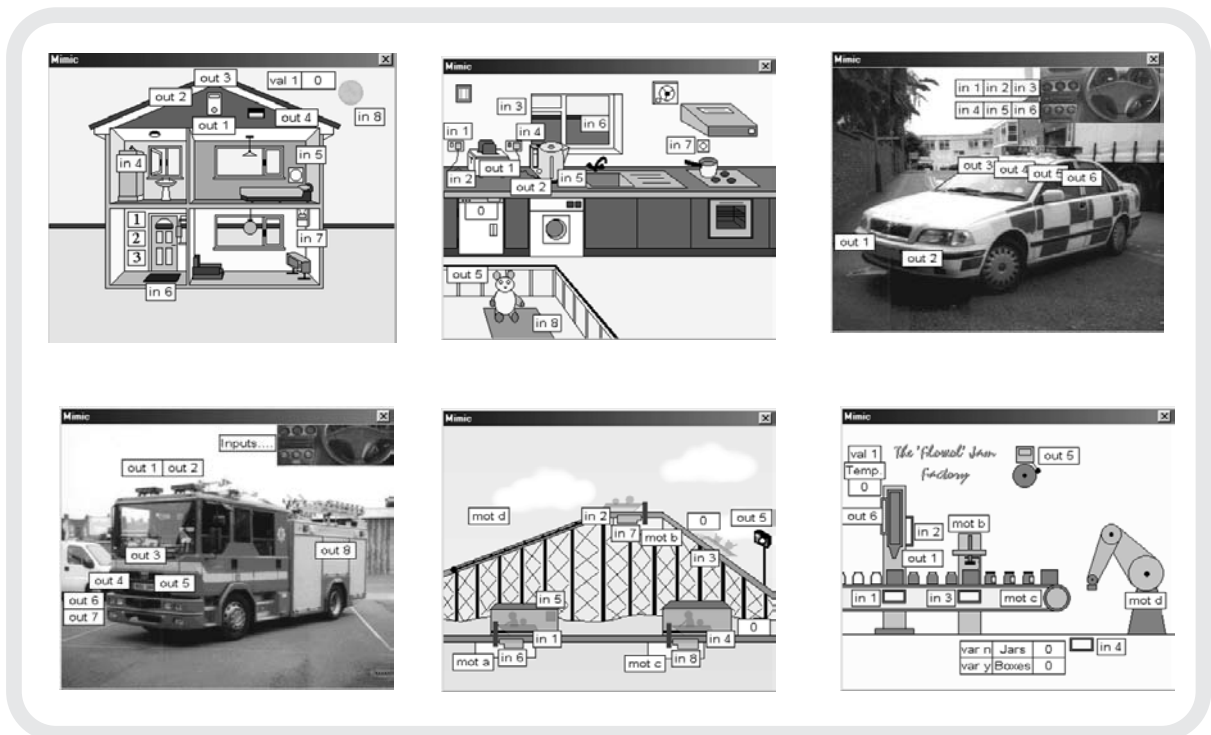

Keying into the National Curriculum

The following tables show how Flowol, FlowGo and the Starter models support both the ICT and Design & Technology Schemes of Work for Key Stages 1 and 2 published by QCA (Qualifications and Curriculum Authority) for England and Wales.



ICT Scheme of Work References

Using Flowol, FlowGo and Model Pack 1 to support the use of the ICT Schemes of Work for Key Stages 1 and 2 published by the Qualifications and Curriculum Authority (QCA)

ICT Scheme of Work

ICT Unit 4E Modelling effects on screen	
What the unit is about :-	How Flowol supports.
Pupils enter instructions to control a screen turtle and compare the operation to a floor turtle. They learn how to write a procedure that 'teaches' the computer a new word and will be asked to produce particular shapes on screen.	Write a program to control one of the Model Pack mimics with Flowol in Simulate mode.
ICT Unit 5E Controlling devices	
What the unit is about :-	How Flowol, FlowGo and the models support.
Pupils learn how to control simple devices, such as buzzers, small motors and lights, using basic control boxes. They learn how to control devices by turning them on and off according to a set of instructions. They will apply this to examples such as traffic lights.	The models in Model Pack 1 include simple devices such as lights and buzzers, which can be turned on and off using a sequence of instructions from Flowol when connected to a FlowGo interface.
ICT Unit 6C Control and monitoring – What happens when	
What the unit is about :-	How Flowol, FlowGo and the models support.
Pupils learn how to use input devices or switches to control a process. They learn that it is possible to attach devices such as pressure pads, light sensors, on/off and magnetic switches etc to a control box. Programs can then be written to carry out a process.	The Lighthouse model has a light sensitive switch. A program can be written so that the lighthouse only flashes when it is dark. The Car Alarm model has three input devices, one on the car door, and 2 push button switches on the key. A program can be written that makes the alarm sound if the car door is opened. The buttons on the key can be programmed to turn the alarm system off and on.
ICT Year 7 Unit 6 Control: input, process and output	
What the unit is about :-	How Flowol, and FlowGo supports.
Pupils learn about control through modelling the working of a carpark barrier. Pupils program a simple cause and effect model e.g. pressing an input switch that produces an output, and then develop their knowledge by using procedures as building blocks.	Using the Car Park Barrier from the Flowol Secondary Mimics pupils can write a program with sub-routines to control the barrier with Flowol in Simulation mode. A 3D model barrier could be built with switches and a motor to be controlled using FlowGo. See also the Fischertechnik Control Kit models available from Data Harvest.

Design and Technology Scheme of Work References



Using Flowol, FlowGo and Model Pack 1 to support the use of the D&T Schemes of Work for Key Stages 1 and 2 published by the Qualifications and Curriculum Authority. (QCA)

Design and Technology Scheme of Work.

D&T Unit 4D Alarms	
What the unit is about :-	How Flowol, FlowGo and the models support.
Pupils develop an understanding of simple electrical control through the designing and making of an alarm system to protect a valuable artefact, which has been brought into the classroom. They are also introduced to the idea of feedback.	An alarm system can be designed for the Car model to protect it from theft, and a Flowol program written to register whether the car door is open or shut by using a proximity sensor as feedback.
D&T Unit 4E Lighting it up	
What the unit is about :-	How Flowol, FlowGo and the models support.
Pupils apply knowledge about electric circuits that they have acquired in science to design and make something that will shine light for a specific purpose. ICT can be used for research and to control their light, switching it on and off when required.	Pupils can create a program for the Fair model, to switch lights on and off to make an attractive display. Pupils can create a program to control the Lighthouse model, using the light switch to detect when it is dark so that the light flashes as a warning to ships at night.
D&T Unit 6C Fairgrounds	
What the unit is about :-	How Flowol and FlowGo supports.
The focus of the unit is for pupils to design and make a model of a fairground ride, but it could be adapted to suit any product in which an electric motor produces rotating movement, which is then transferred using belts and pulleys.	Motors, when connected to FlowGo, can be controlled both forwards and backwards at various speeds. This enables pupils to design a ride that starts slowly, and builds up to full speed, before slowing down again.
D&T Unit 6D Controllable Vehicles	
What the unit is about :-	How Flowol, and FlowGo supports.
As an alternative to fairgrounds, the focus is for pupils to design and make a model of a vehicle in which an electric motor produces rotating movement, which is then transferred using belts and pulleys, then can be controlled using computer control.	Up to two motors, when connected to FlowGo, can be controlled both forwards and backwards at various speeds.



Key Skills

Where Flowol and FlowGo support the development of the 6 Key Skills.

Communication:

- When working with a partner or in small groups, pupils will need to communicate their ideas verbally, in writing, and then in the form of an on-screen flowchart. Pupils could then explain (verbally) what they have done, and how they have approached a problem by talking to a larger group or the whole class.
- Pupils will develop and refine ideas and make things happen by planning and giving instructions to make a device work (*e.g. when devising a procedure to operate lights*).

Application of number:

- Pupils will work with mathematical models when producing flowcharts and sequences of instructions for controlling mimics and actual components/models.
- Pupils will consider and use units of time (seconds and parts of seconds), for example when planning how long devices should stay 'on' or 'off'.

Information technology:

- Through the programmes of study for ICT pupils will have opportunities to develop their knowledge, skills and understanding in relation to developing ideas and making things happen. Specifically pupils will have opportunities
- 'To create, test, improve and refine sequences of instructions to make things happen and to monitor events and respond to them.'
- 'To use simulations and explore models in order to answer 'What if ... ?' questions, to investigate and evaluate the effect of changing values and to identify patterns and relationships.'

National Curriculum (2000)

Working with others:

- When working in pairs or as part of larger groups, pupils will have opportunities to develop their social skills: specifically cooperation, and consideration of the views of other members of the group.

Improving own learning and performance:

- When planning and producing sequences of instructions for controlling mimics or actual components, pupils will review, modify and evaluate their work as it progresses
- Through systematic use of the program, pupils will develop their knowledge, skills and understanding in many aspects of ICT, and also in D&T when actual working products are controlled.

Problem solving:

- Pupils will have opportunities to consider and work through realistic situations and develop a variety of solutions to problems (*e.g. devising a car alarm which triggers when the car door is opened*). This includes: developing skills to identify and understand a problem (*e.g. because of the risk of theft, cars need to be effectively protected*); planning ways to solve the problem (*different options inc - pressure pads, light sensors*); monitoring progress in tackling the problem (*is my solution working yet?*); and reviewing solutions to the problem (*how well does my solution work?*).

The Fair Model and FlowGo






Teacher's Notes for the Fair Mimic & Model

Resources:

- Flowol software
- Fair mimic from the Model 1 Mimic disc
- Fair model
- FlowGo control interface

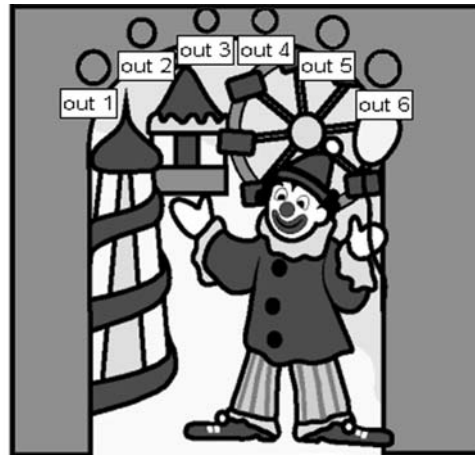

Software checklist:

1. Load Flowol.
2. Select the *Mimic* option from the *Window* menu. Select the **p1-fair** mimic and ensure the **Show Labels** box is ticked.
3. Ensure that the selected interface is FlowGo (using *Interface* from the *Control* menu).
4. If you are using the Fair mimic only, stay in **Simulation** mode. If you are using with the fair model attached to FlowGo, untick **Simulate**.

Fair Model

The Fair model has a display of 6 lights above the entrance to a fairground. Connect the model to FlowGo and check that the numbering of the outputs matches those on the mimic.

Learning objectives/outcomes

Pupils will learn:

- That a computer can be given instructions (programmed) to carry out a sequence of actions.
- How to create a program that turns outputs on and off.
- How to create a program that runs continuously.
- That a program can be edited.

Extension activities – Pupils will learn about the use of sub-routines.

Activities

Activity 1 Explore the Fair mimic





- Show how the outputs respond on the mimic and the model.

Activity 2 Switching one light on and off.

- Introduce the Flowol symbols and what they do.
- Show how to select and place the symbols on the screen.
- Show how to select the instructions from the Prompt Box.
- Show how to link the symbols together using the Arrow tool.
- Show how to run a program.
- Show how to save this program, and call it Flash.



Teacher's Notes for the Fair Mimic & Model



To select a symbol from the toolbar		Use the mouse to click once on the symbol (click and release). Move the mouse to the right. The symbol will appear in the workspace, click to place it in the required position.
	Start	This symbol is used to start or stop a program.
	Output or Action symbol	This symbol is used to switch outputs on or off.
	Process (Delay) symbol	This symbol is used to make the program wait (delay) before going on to the next instruction.
	Arrow or Line	This arrow is used to link the symbols together to show the order in which each symbol (instruction) should be carried out

Activity 3 To make the light flash all the time

- Erase the Stop symbol. Show how to use the 'Helping Hand' to highlight a symbol, and how to erase.
- Click on the arrow (the line tool) and connect the Delay 2 at the bottom of the flowchart to the first action symbol (Turn Output 1 on)
- Run the program and ask the pupils to describe what is happening.
- Emphasise that the program is running in a LOOP.
- Save the changed program as *Flashes*

	The Helping Hand	Click on the helping hand, click on the item in the flowchart (it will change colour).
	Erase tool	After using the helping hand to select the item, click on erase to delete.

Activity 4 To turn more than one light on and off at the same time.

- Show how to use the 'Helping Hand' tool to highlight the symbol to be changed
- (Turn Output 1 on).
- In the Prompt Box, click on the Clear button, and put in new instructions (Turn Output 1 on, 6 on, then OK).





Teacher's Notes for the Fair Mimic & Model

Challenges

A sample program is included for each challenge on the Activity Sheet. These samples are for guidance only, as there is always more than one way to write the program!

Challenge 1: Change the delay instructions so that the light will stay on longer than it is off.

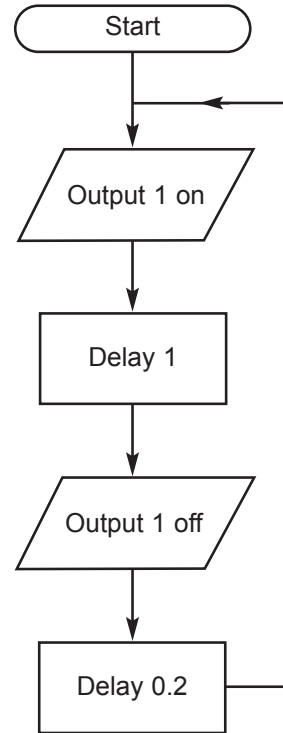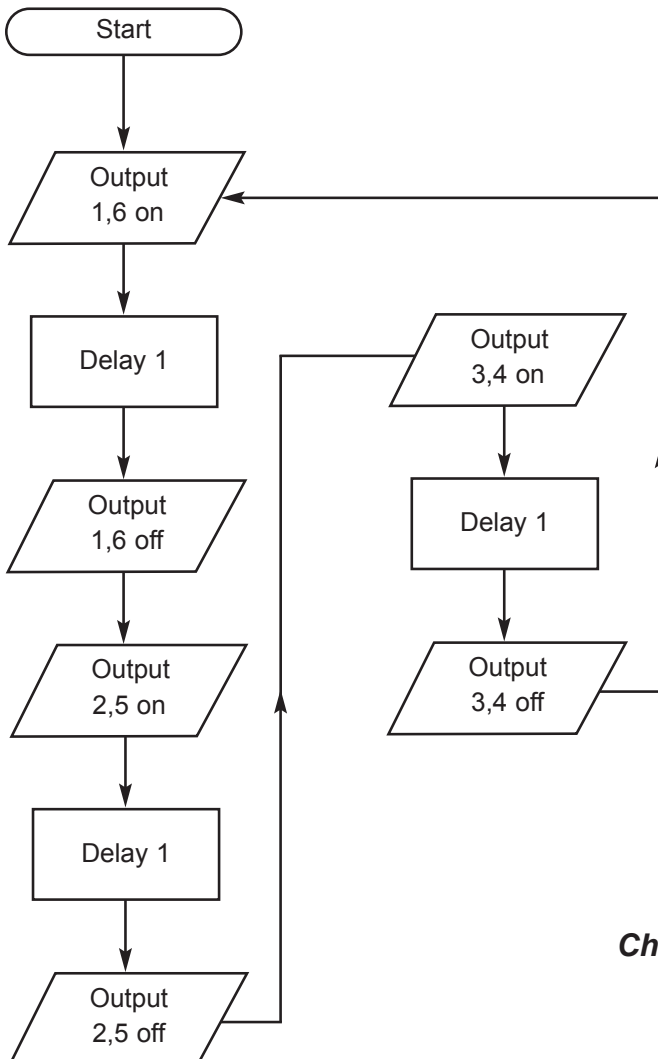

Challenge 2: Create a program to flash the two red lights on and off, then the two green lights on and off, and then the two yellow lights on and off.

Challenge 3: Create your own light display.

Tip: A Delay can be as short as 0.1 second.



Teacher's Notes for the Fair Mimic & Model

Extension Activity Challenges:


Return to the Fair model after introducing sub-routines with the Lighthouse model for the final challenges.

1. Using sub-routines to repeat a set of instructions

Instead of writing a long list of on/off instructions to flash the lights, a sub-routine can be created to flash the light once, called *Flash1*. This *Flash1* sub-routine can be then be called from the main program and run 3 times.


Step 1: Create a sub-routine

- Select the Start Symbol
- In the Prompt Box, click on Sub.
- Type the name *Flash1* in the text box, and then click on OK.
- Copy the instructions shown

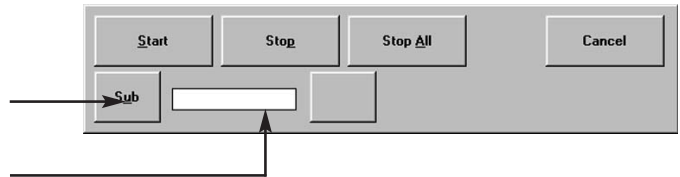

Step 2: Using the sub-routine in the main program

- Select Start
- Select the Process (Delay) Symbol.
- Click on Sub and select 1(Flash1)
- Click on the number 3, OK
- Select Stop
- Link up the flowchart.

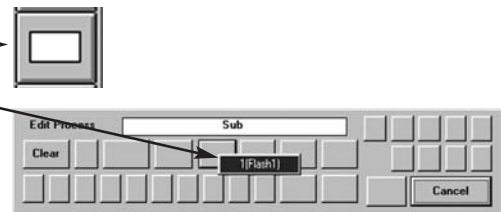


Save this program as *3flashes*

Ask the pupil to run the program and watch the mimic.

Ask if the lights flash 3 times, then ask them to watch the symbols and observe how the program performs.

Challenge 4: Create a flashing light display with the red lights flashing twice, then the yellow lights 3 times, then the green lights twice - make it run continuously.

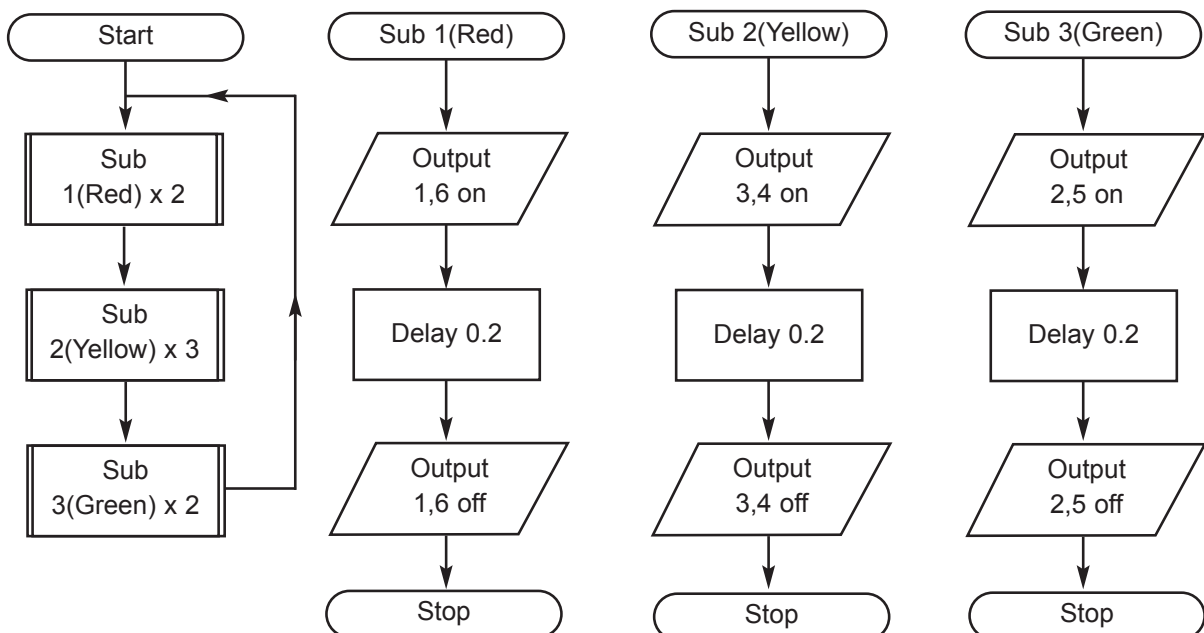



Teacher's Notes for the Fair Mimic & Model

Challenge 5: Program the lights to represent a ball being thrown from the Clown's left hand to right hand, and then back again.



Teacher's Notes for the Fair Mimic & Model



2. Simple on/off activities using other available mimics (see page 33)

Bridge Lights - Create a traffic light sequence.

Robot Face - Make the eyes open and shut. Make the mouth open and shut, and make the aerial rise and lower.

School Patrol Crossing - Make the sign flash.

Police Car and Fire Engine - Make the lights flash.

Clown Face - Make the eyes wink and then blink.

Traffic Lights - Create a traffic light sequence.

Cross-curricular links:

Literacy - Stories about fairgrounds

Science - Simple electrical circuits

Design & Technology (Primary) - Designing any display of lights e.g. Celebrations, Christmas tree lights, Festivals, Diwali, Bonfire night, Advertising displays, Street illuminations, creating a book display, etc.



Activity Sheet

Activity 1:


Explore the Fair mimic
Click on the lights over the fairground entrance to turn them on and off.



Activity 2:

Switching one light on and off
You are going to create a simple program to switch the light on, keep it on for 2 seconds, and then switch it off.


To run your program, click on the red RUN button.
Save your program and call it *flash*.

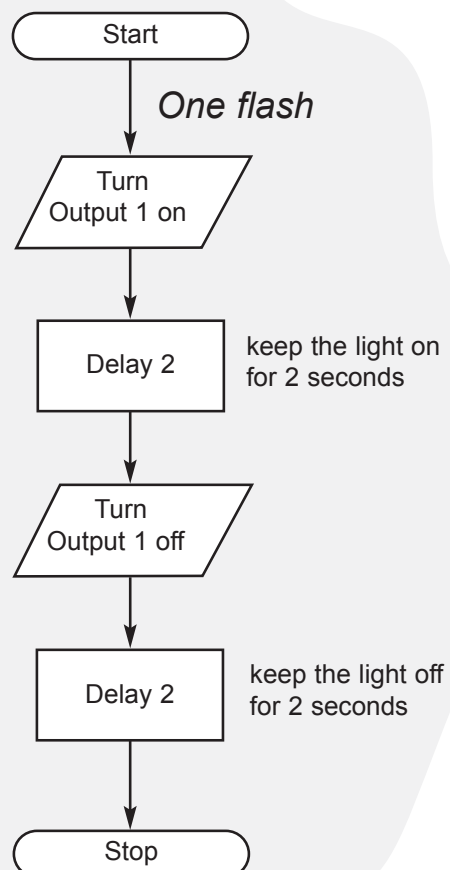



Activity 3:

To make the light flash all the time
Change your program *flash*. See below: -

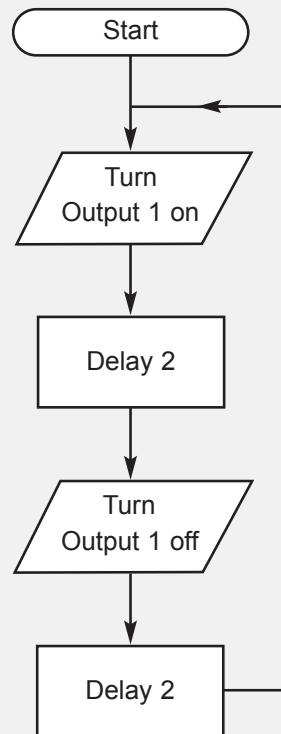

Run your program. What happens?

You will have to click on the red STOP button to stop the program running.

Save the changed program as *flashes*.

Challenge 1:


Change the delay instructions so that the light stays on for longer than it is off.



Activity Sheet

Activity 4:

To turn more than one output on or off at the same time
Change your program so that it looks like this: -

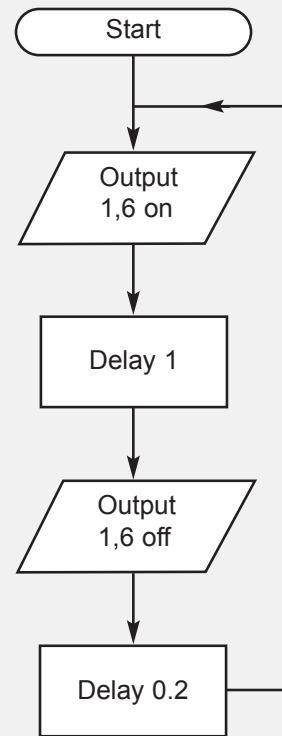

Challenge 2:

Can you create a program to flash the two red lights on and off, then the two green lights on and off, and then the two yellow lights on and off? It should run all the time.

Challenge 3:

Create your own light display.

Extension Activity Challenges:

Challenge 4:

Create a flashing light display that flashes the red lights twice, the yellow lights 3 times, then the green lights twice. The program should run all the time.

Challenge 5:

Can you program the lights to represent a ball being thrown from the Clown's left hand to right hand and then back again?

The Lighthouse Model and FlowGo






Teacher's Notes for the Lighthouse Mimic & Model

Resources:

Flowol software
Lighthouse mimic from the Model 1 Mimic disc
Lighthouse model
FlowGo control interface

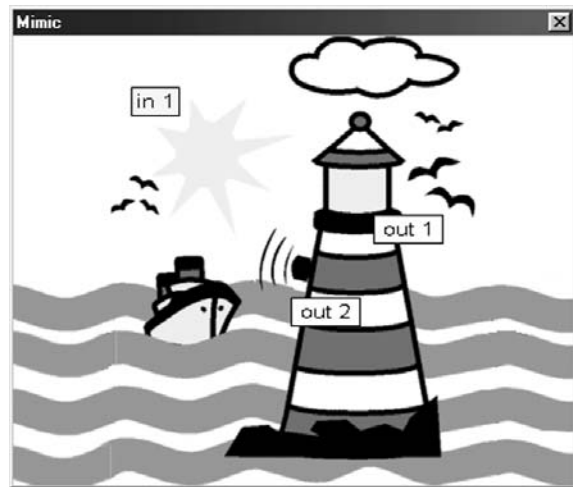

Software checklist

1. Load Flowol.
2. Select the *Mimic* option from the *Window* menu. Select the **p1-light** mimic and ensure the **Show Labels** box is ticked.
3. Ensure that the selected interface is FlowGo (using *Interface* from the *Control* menu).
4. If you are using the Lighthouse mimic only, stay in **Simulation** mode. If you are using with the Lighthouse model attached to FlowGo, untick **Simulate**.

Lighthouse Model

On the Lighthouse model, the lighthouse lamp is Output 1 and Output 2 is the foghorn (buzzer). There is a light sensitive switch (Input 1) that is used to detect whether the sun is shining or not. Connect the model to FlowGo and check that the numbering of the outputs/inputs matches those on the mimic.

Learning objectives/outcomes

Pupils will learn:

- that a computer can be given instructions (programmed) to carry out a sequence of actions.
- that feedback from an input signal can be used in a program to decide whether to carry out a set of instructions or not.
- that a program can be divided up into sub-routines and repeated for a set number of times.
- to develop a system that can control events in response to conditions.

Activities



Activity 1 Exploring the lighthouse model and mimic

- Show that clicking on the mimic switches on and off outputs to both the mimic and the model.
- The yellow indicator LED situated next to Input 1 on FlowGo will be lit when light is reaching the light sensitive switch. Show the change to the input when the light switch is covered up i.e. when it becomes dark.

Teacher's Notes for the Lighthouse Mimic & Model




Activity 2 Using a switch as an input signal


	Decision box	Used to check for feedback from an input signal
	Yes / No Lines	Where a decision box has been used, then both a Yes and a No line have to be selected.

Introduce the Decision symbol and demonstrate the Yes and No logic. The example program is annotated to help explain the decision process.

- If the sun is shining (the input is ON), then the program keeps looking again and again at the state of the input.
- If the sun is not shining (the input is OFF), then the program follows the instructions to flash the lighthouse lamp.

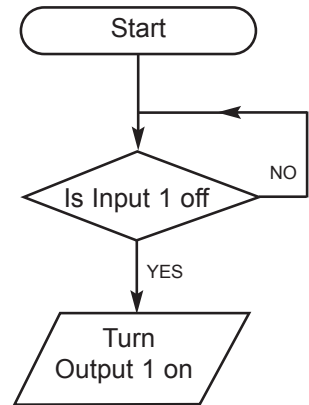

See the following examples and the Flowol Tutorial for more information on using the Decision symbol.

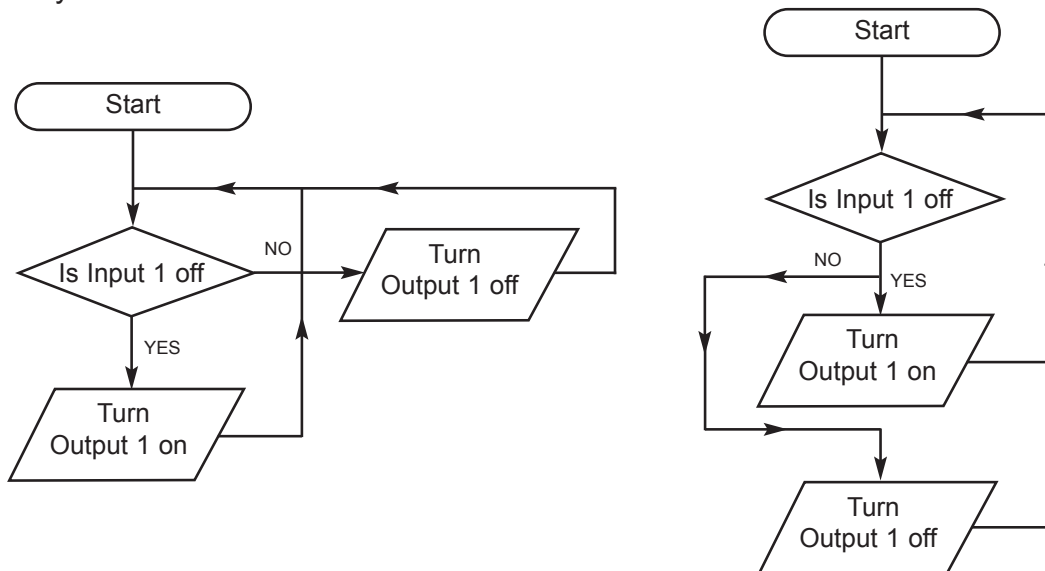

Activity 3 Making the lighthouse repeat a flashing light pattern - introducing sub-routines.

The purpose of this activity is to show pupils how to break up a very long string of instructions into short, manageable sub-routines. Sub-routines can be called from the main program and instructed to run just once or repeated a set number of times.




Teacher's Notes for the Lighthouse Mimic & Model

The example programs the lighthouse lamp to give 2 long flashes then 3 short flashes, whenever it is dark.

Step 1: Create a sub-routine

- Select the Start Symbol
- In the Prompt Box, click on Sub.
- Type the name *long* in the text box, and then click on OK.

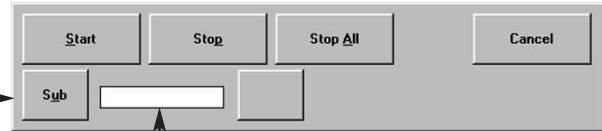


Step 2: Create a second sub-routine

- Select the Start Symbol
- In the Prompt Box, click on Sub.
- Type the name *short* in the text box, and then click on OK.

Step 3: Create the main program

When selecting the sub-routines:

- Select the Process (Delay) Symbol
- Select Sub 1(long) then click on the number 2 and OK.
- Select Sub 2 (short) then click on the number 3 and OK.
- Link up the flowchart as shown.

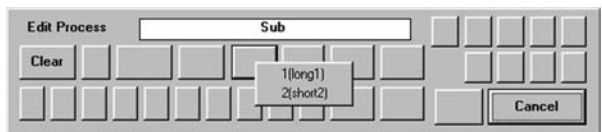

Save this program as *Rocky*.

Ask the pupils if the lighthouse flashing sequence looks effective. Then ask them to watch the symbols and observe how the program performs.

Teacher's Notes for the Lighthouse Mimic & Model




Extension activities/Challenges

A sample program is included for each challenge. These are intended as guidance only, there is always more than one way to write a program!


Challenge 1 :

Create a simple program that flashes the light on and off all the time. This program will utilise the skills learned in the Fair Activity.

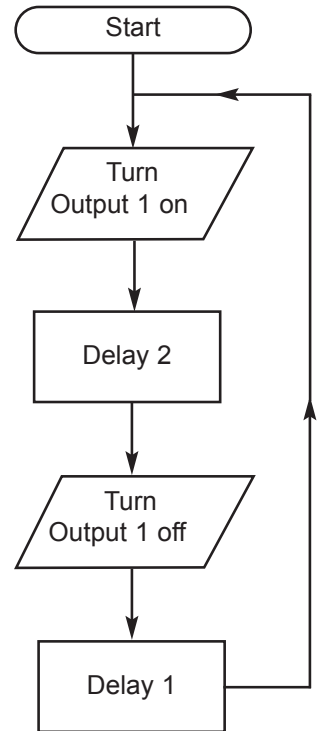

Challenge 2:

Create a new program to flash the lighthouse when the sun is not shining. The light sequence should be 1 long flash then 2 short flash.

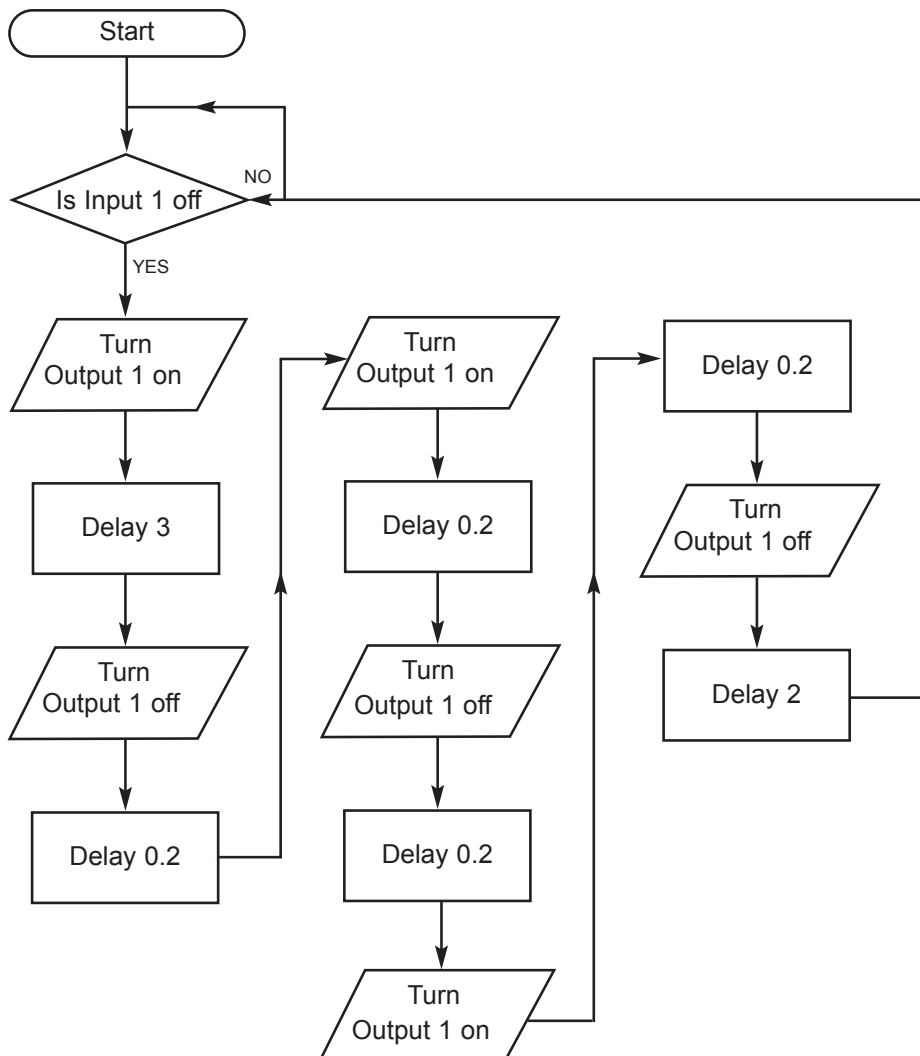




Teacher's Notes for the Lighthouse Mimic & Model

Challenge 3:

Create a program that, when the sun is not shining, sounds the foghorn once after every 5 flashes of the lighthouse.

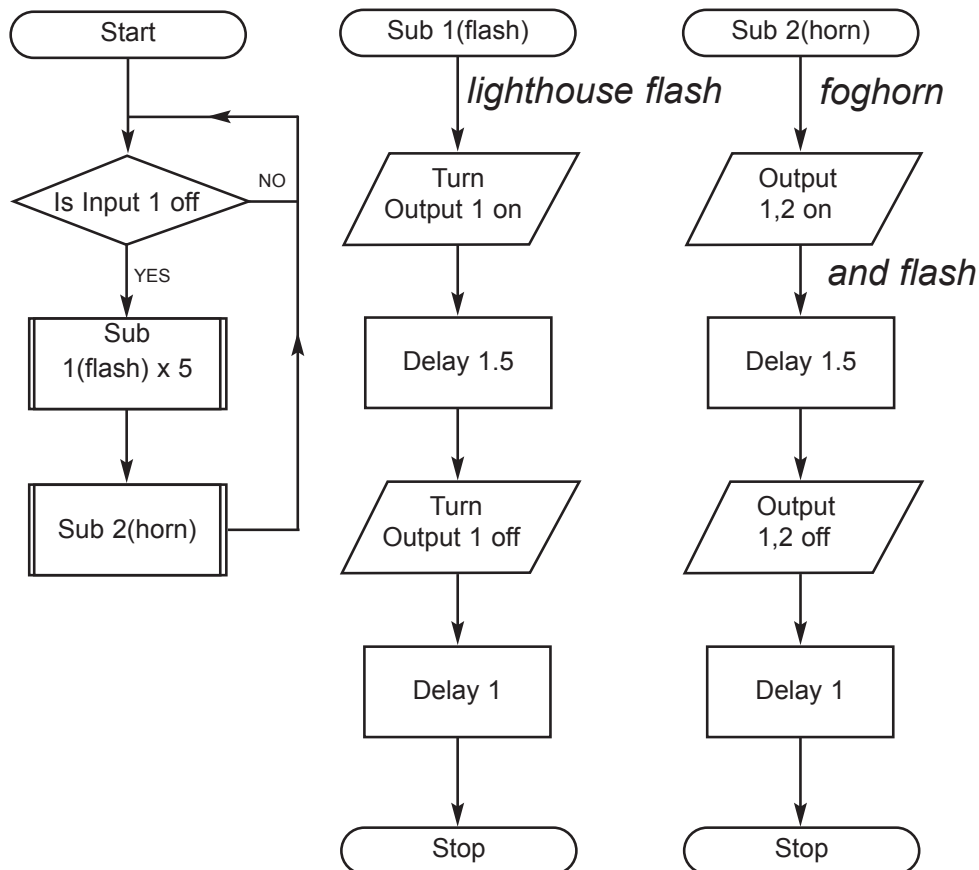

Other Extension Activities:

Using sub-routines in programs using other available mimics (see page 33).

Pelican Crossing – flashing amber light

Robot Face - making the eyes wink and the mouth open and shut

Simple Level Crossing – flashing lights

Fire Engine – flashing lights

Simple Home Alarm System – flashing light as part of the alarm sequence

Cross-curricular links:

Literacy - Lighthouse Keeper stories

Science - Making simple circuits with a bulb and a switch.

Night and day studies


Design & Technology (Primary) - Designing with found or reclaimed materials and incorporating a simple electrical circuit e.g. vehicle with warning light, torch/beacon with a light sensitive switch, advertising display that is triggered by a pressure mat etc.



Activity 1:

Exploring the lighthouse model and mimic


Try clicking on the lighthouse lamp and beside the foghorn on the mimic. What happens? Cover the light sensitive switch (sun) on the model. Does the LED next to Input 1 on FlowGo change when the switch is covered?

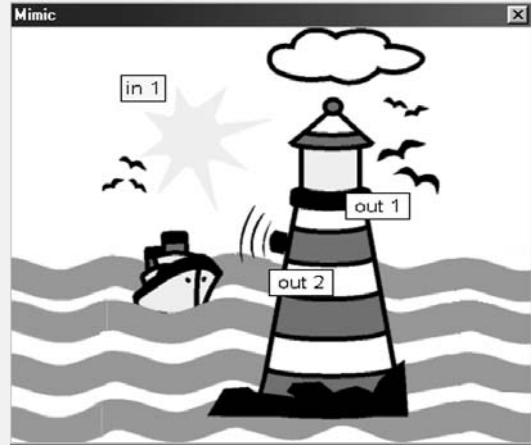

Challenge 1:

Create a simple program that flashes the light on and off all the time.

Activity 2: Using a switch as an input signal.

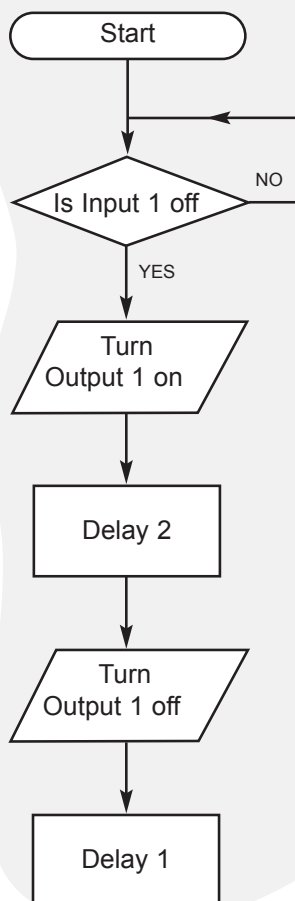

Is Input 1 Off
(Is it dark? Is the sun NOT shining?)

No - Its still light (the sun is shining) so keep checking.

Yes - run the instructions to flash the lamp.

Save this program and call it *Dark1*.



Activity Sheet

Challenge 2:


Change your program to flash the lighthouse when the sun is not shining. The light sequence should be 1 very long flash then 2 short flashes.

Activity 3 Making the lighthouse repeat a flashing light pattern - introducing sub-routines.

Create sub-routine long1 first.

Then create sub-routine short2.

The main program should run sub-routine 1 once, and repeat sub-routine 2 twice.

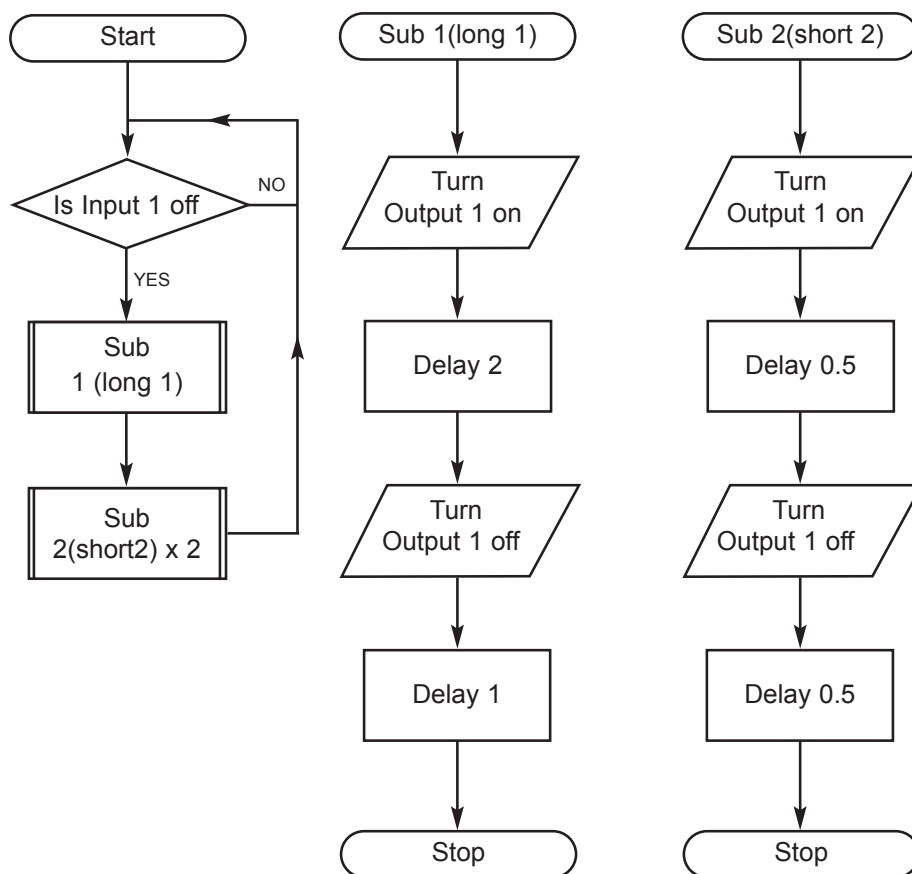


Challenge 3:

Can you create a program that sounds the foghorn once after every 5 flashes of the light, only when the sun is not shining?

The Car Model and FlowGo




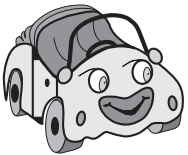

Teacher's Page for the Car Alarm Mimic & Model

Resources:

- Flowol software
- Car mimic from the Model 1 Mimic disc
- Car model
- FlowGo control interface

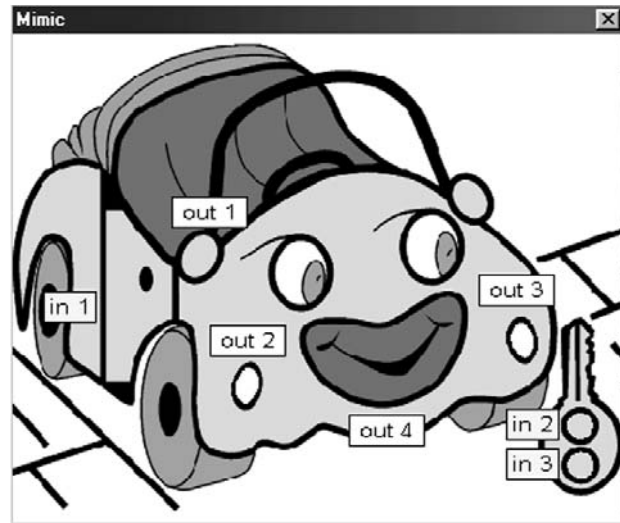

Software checklist:

1. Load Flowol.
2. Select the *Mimic* option from the *Window* menu. Select the **p1-Car** mimic and ensure the **Show Labels** box is ticked.
3. Ensure that the selected interface is FlowGo (using *Interface* from the *Control* menu).
4. If you are using the Car mimic only, stay in **Simulation** mode. If you are using with the car model attached to FlowGo, untick **Simulate**.

Car Model

The Car model has:

- Light on wing mirror – Output 1
- Left Headlight – Output 2
- Right Headlight – Output 3
- Horn – Output 4
- Car Door switch – Input 1
- Key button – top – Input 2
- Key button – bottom – Input 3

Learning objectives/outcomes

Pupils will learn:



- that a computer can be programmed to monitor a system.
- how to use an input signal in an alarm system.
- that a complex system can be broken down into manageable parts.
- to develop a system that can control events in response to conditions.
- to use input and output devices and produce a simple set of instructions linking causes and effects.

Activities

Activity 1 Exploring the car model and mimic

- Show that clicking on the mimic causes outputs to switch on and off on both the mimic and the model.
- Show the effect of opening and shutting the car door (Input 1 is on when the door is shut so the yellow LED on FlowGo will be lit).

Teacher's Page for the Car Alarm Mimic & Model

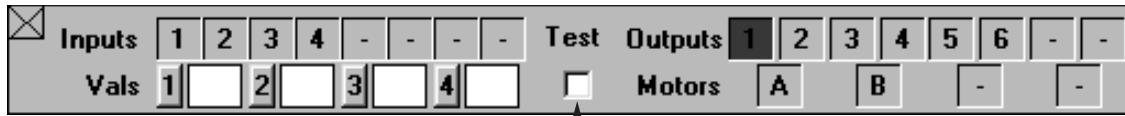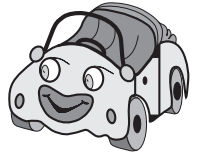

Test box


To show the change in state of the inputs (with the model connected to FlowGo):


- Open the **Small Window** (from the Window menu). Click so a tick appears in the **Test box**.
- Press the button of the push-to-make switches on the key (to switch it ON) and watch the box below the relevant input light up.
- The box below Input 1 (the car door) will be lit if the car door is shut (Input On).

Activity 2 Using more than one flowchart

This activity introduces the concept of an emergency stop. Using the Stop All command in the second flowchart stops the program running and switches everything off.

Flash the lights and sounds the horn if the car door is opened

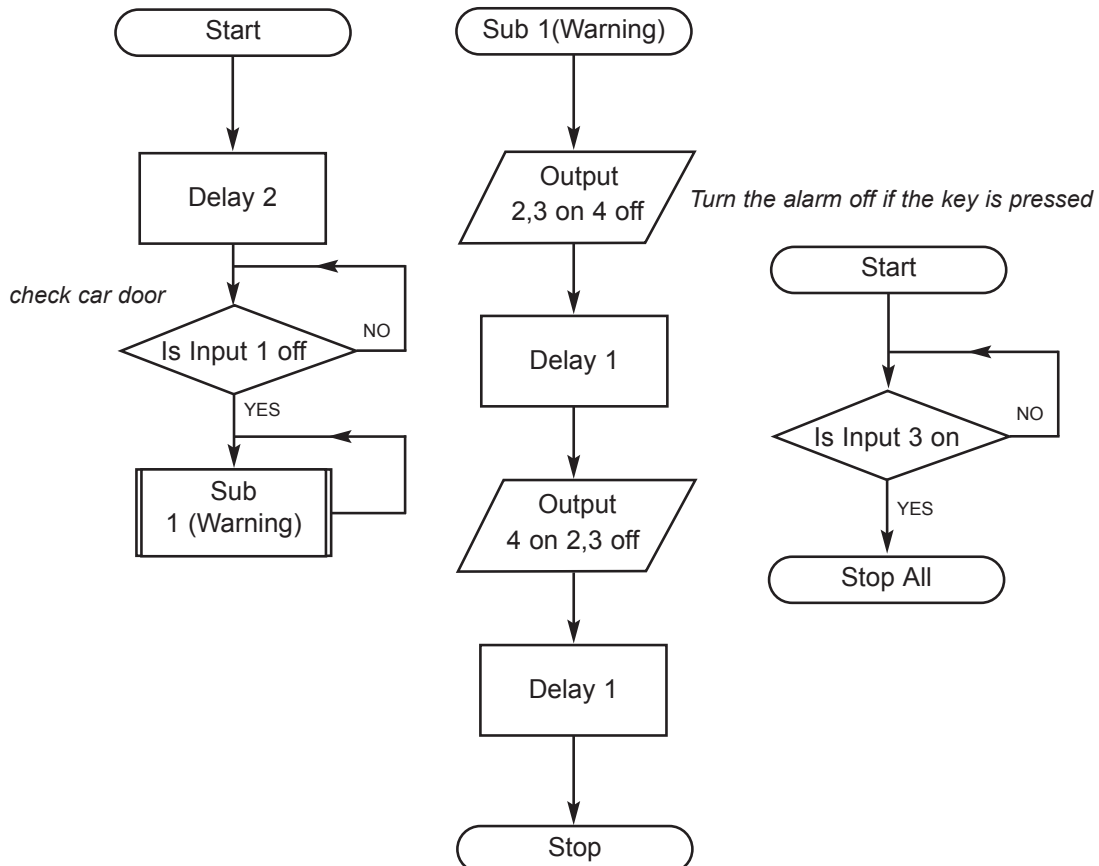

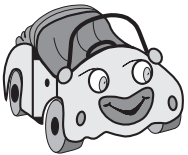

Teacher's Page for the Car Alarm Mimic & Model

Activity 3 Planning a complete car alarm system.

This activity guides pupils through breaking down a control system into separate actions (sub-routines). This program should be planned out first on paper. This plan can then become the design specification, against which the final program can be appraised.

Example

Plan of the Car Alarm system

- The program should continuously check to see if an alarm has been set (has the top switch on the key been pressed?).
- When the alarm has been set pressing on the bottom switch on the key should deactivate it.
- The alarm needs to activate if the door is opened.
- If the alarm has been triggered it should continue until the bottom switch on the key has been pressed.
- Shutting the car door should not deactivate the alarm.

Sub-routines needed (these should be created first):


	Sub-routine name
Flash the indicators once (turn on, one second delay, turn off).	Flash
Sound the buzzer and flash the indicators to show the alarm has been triggered (on, half a second delay, off).	Alarm
If the door is open activate the alarm (ALARM). If the door is closed, indicate that the alarm has been set by the indicators flashing twice (FLASH x 2).	Set It
Turns the alarm status light on. Checks the door. <ul style="list-style-type: none">• If the door is shut check for the bottom switch on the key to be pressed (DISARM).• If the door has been opened activate the alarm (ALARM). Keep the alarm running even if the door is then shut until the bottom key is pressed (DISARM).	Armed
Show the alarm is switched off by one long flash of the indicators (FLASH) and switching the alarm status light off.	Disarm

The sub-routine Armed may need some explanation, as there are two inputs to check.

1. If the door is open (Input 1 Off) then run the Alarm sub-routine, until the Alarm Off button (Is Input 3 On) is pressed.
2. If the Alarm Off button is pressed (Is Input 3 On) then disarm the system.

The main program should continuously check to see if the top switch on the key (input 2) has been pressed and then activate the alarm.


Teacher's Page for the Car Alarm Mimic & Model

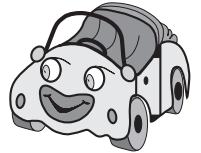

Extension Activities/Challenges

Challenge 1:

Create a program that flashes both indicator lights on and off. This is a simple activity using the same skills as learned in the Fair activity.

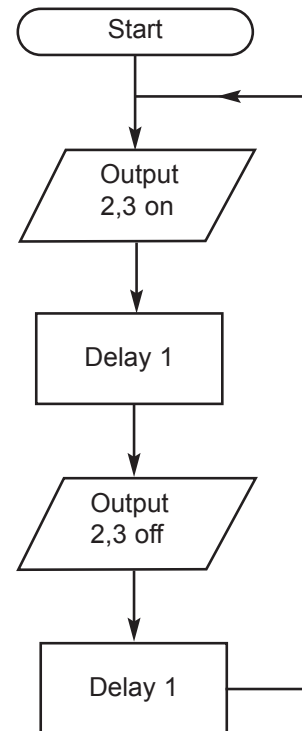



Challenge 2:

Create a simple alarm program. The indicator lights should flash and the horn sound if the car door is opened. The alarm should continue, even if the door is closed again.

Flash the lights and sound the horn if the car door is opened


The solution should include the use of a sub-routine (as in Activity 3 of the Lighthouse) but there is another solution that works - see Car Challenge 2b

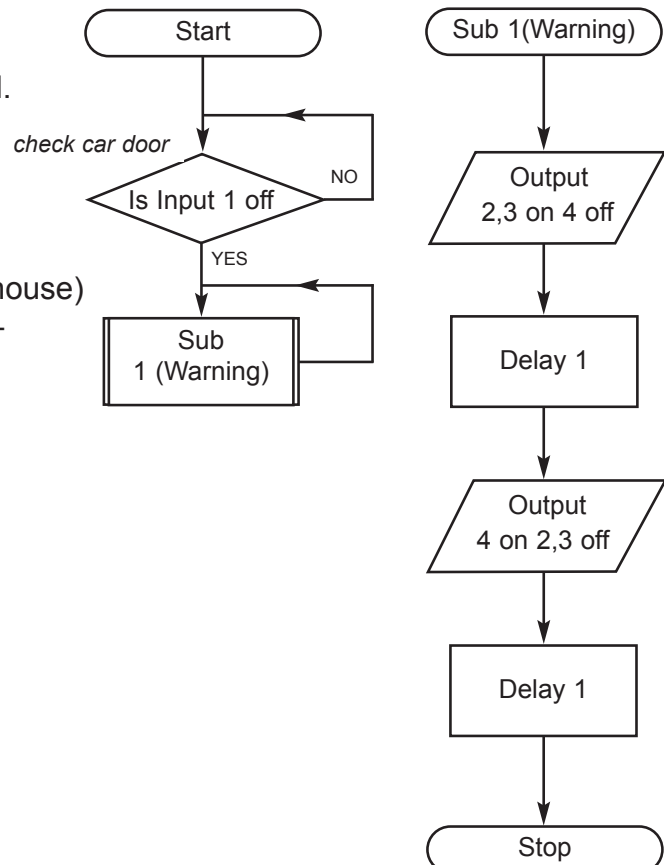

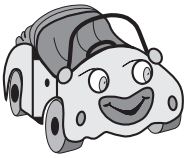

Teacher's Page for the Car Alarm Mimic & Model

If pupils produce this type of program, show them how to edit their program to create a sub-routine for the flashing lights and horn (Alarm). Pupils will need to be able to use a sub-routine in order to complete the other activities.

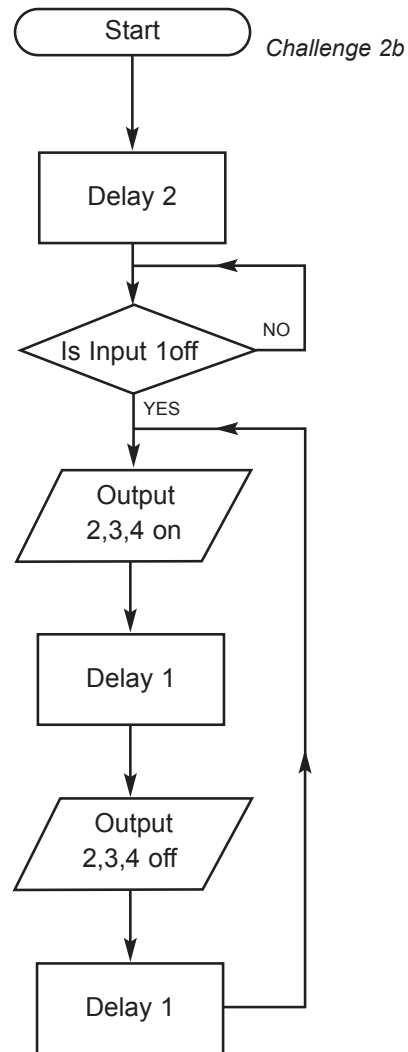

Challenge 3:

Can your Car Alarm program be improved so that it can respond instantly to the Alarm Off button?


Cross-curricular links:

Literacy - Crime stories.

Numeracy - Logic.

Design & Technology - Design a house alarm using a cardboard box and a selection of switches to detect windows and doors being opened, or some one standing on a mat.


Design an alarm system for a valuable item e.g. jewel, museum exhibit, etc.

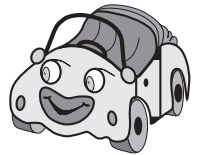

Activity 1:

Exploring the Car model and mimic

Use the mouse to click on the indicator lights (output 2 and 3) and switch them on and off. Click on the wing mirror (output 1) to switch the alarm status light on and off.

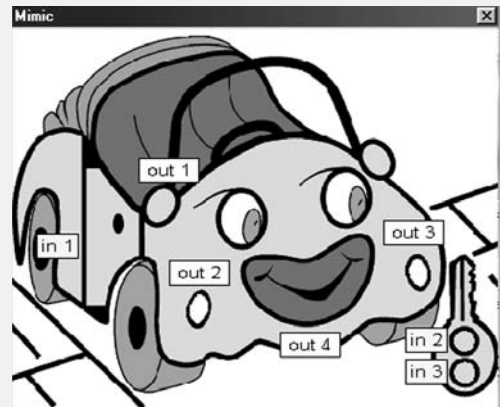

From the Window menu select the Small Monitor and ensure Test box is ticked.

<input checked="" type="checkbox"/>	Inputs	1	2	3	4	-	-	-	-	Test	Outputs	1	2	3	4	5	6	-	-
	Vals	1		2		3		4		<input checked="" type="checkbox"/>	Motors	A	B	-	-	-	-	-	-

Press the push buttons on the key and see the change on the mimic and in the box below input 2 and 3 on the small monitor window.

Open the car door on the model. Look at input 1 on FlowGo and the box below input 1 on the small monitor window. Is the car


Challenge 1:

Create a program that flashes both indicator lights on and off. Save your program as *car1*.

Challenge 2:

Create a simple alarm program. The indicator lights should flash and the horn should sound if the car door is opened. The alarm should continue, if the door is closed again.

Tip! If using the mimic only, remember when you run your program shut the car door before setting the alarm. Save your program as *car2*.

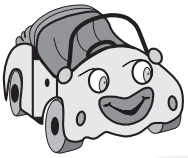

Activity Sheet

Activity 2:

Using 2 flowcharts

The lower button on the key (input 3) is used to turn the car alarm off.

Flowol can run more than one flowchart at the same time. The second flowchart is used to check if the key button (input 3) has been pressed.

The command Stop All is used to stop the whole program running and switch everything off.



Add this second flowchart to your program, *car2*.

Run the program and trigger the alarm by opening the car door.


Press the lower button on the key. Did the alarm stop?


Save this program as *carkey*

Flash the lights and sound the horn if the car door is opened

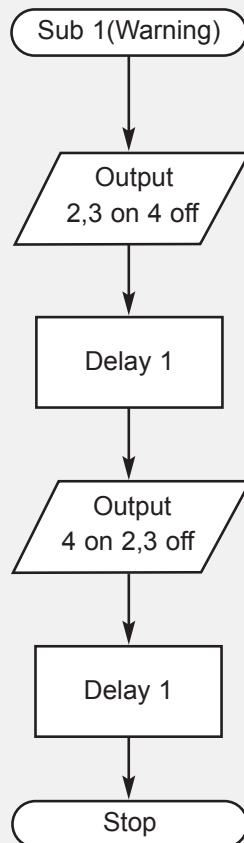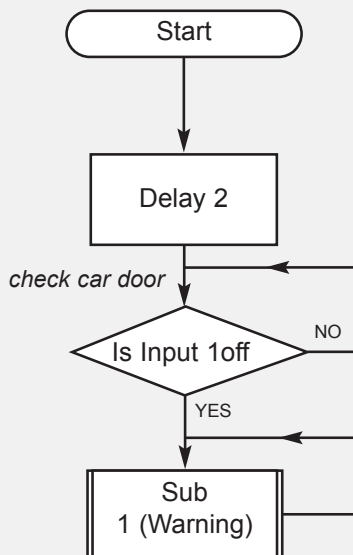

Turn the alarm off if the key is pressed

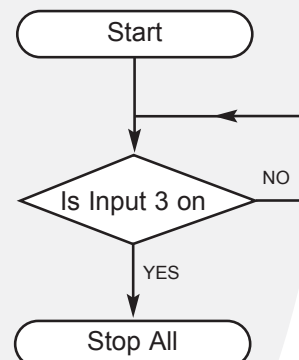

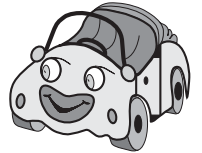

Activity 3:

Planning a complete car alarm system.

The ultimate challenge is to design a program that is running all the time.

This program shows a car alarm system that is set by pressing the Alarm On button on the key, and is disarmed by pressing the Alarm Off button. Once the system is set, opening the car door triggers the alarm, making the horn sound and the indicator lights flash.


The alarm can only be switched off by pressing the Alarm Off button.

To plan this program, write a list of the sequence of actions required.


Example:

	Sub-routine Name	Sub-routine Action
1	Flash	Single flash of both indicators with a 1 second delay to be used when alarm is set, and when the alarm system is switched off.
2	Alarm	Sound the horn and flash both indicators on and off once with half a second delay.
3	Set it	Check if door is shut, run Flash twice, then run Armed.
4	Armed *See example	Turn the alarm status light on. Check the door <ul style="list-style-type: none"> • If opened, run alarm then check key button and if it is pressed, run Disarm. Keep running alarm if key is not pressed. • If the door is not opened, check key for button press, if not pressed, recheck door.
5	Disarm	Run Flash once, and turn the status light off.


Main Program - checks for input 2 to be ON then sets the alarm.

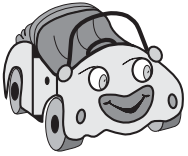

Activity Sheet

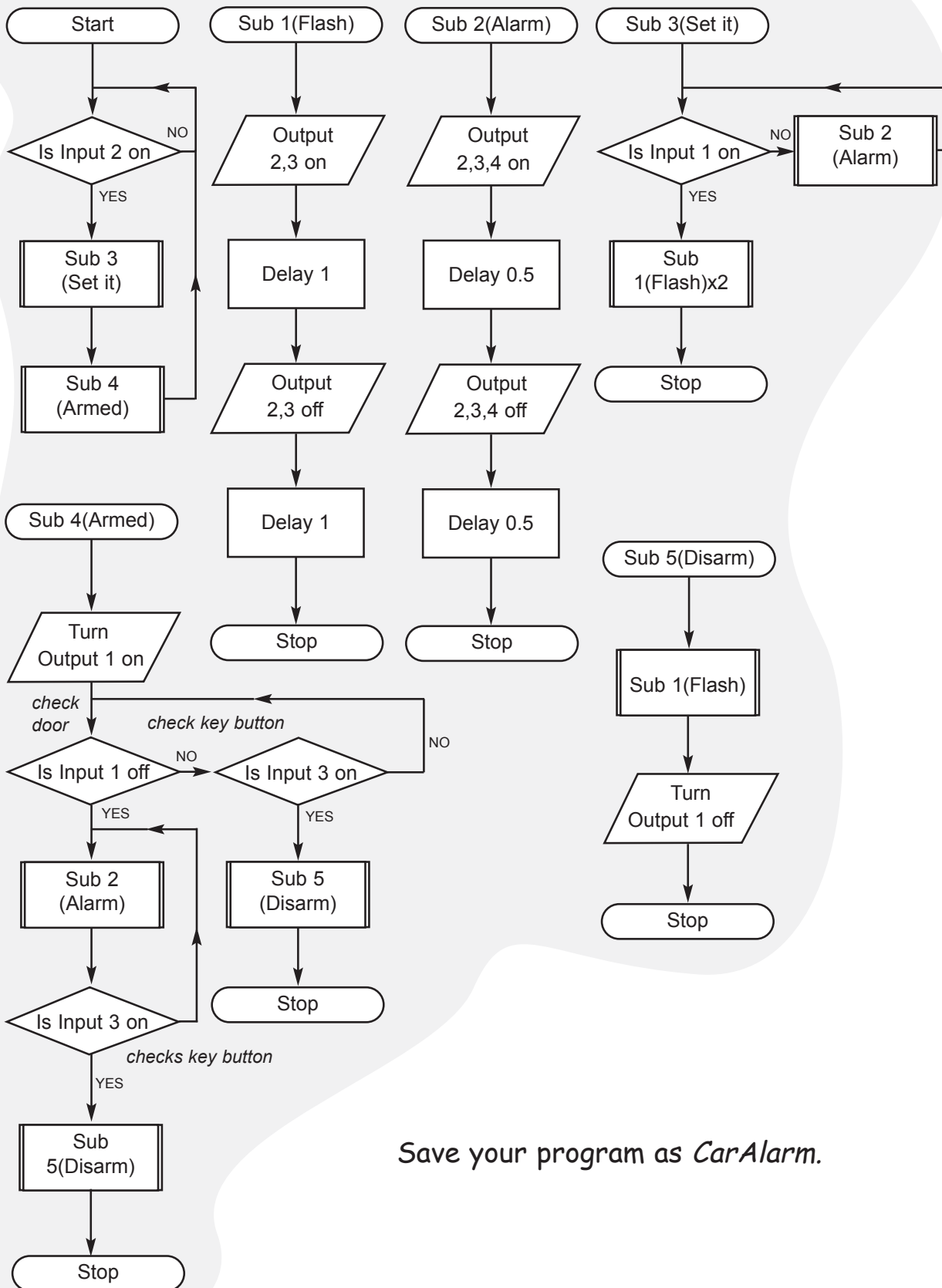

Save your program as *CarAlarm*.

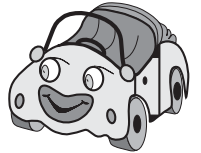

Try your program.

You may find that if the alarm has been triggered you have to hold down the Alarm Off button for some time to stop the alarm system. Can you explain why?

Challenge 3:

Can you improve on your Car Alarm program so that it can respond instantly to the Alarm Off button?